
50X1-HUM

Page Denied

Next 1 Page(s) In Document Denied

MODEL 93B

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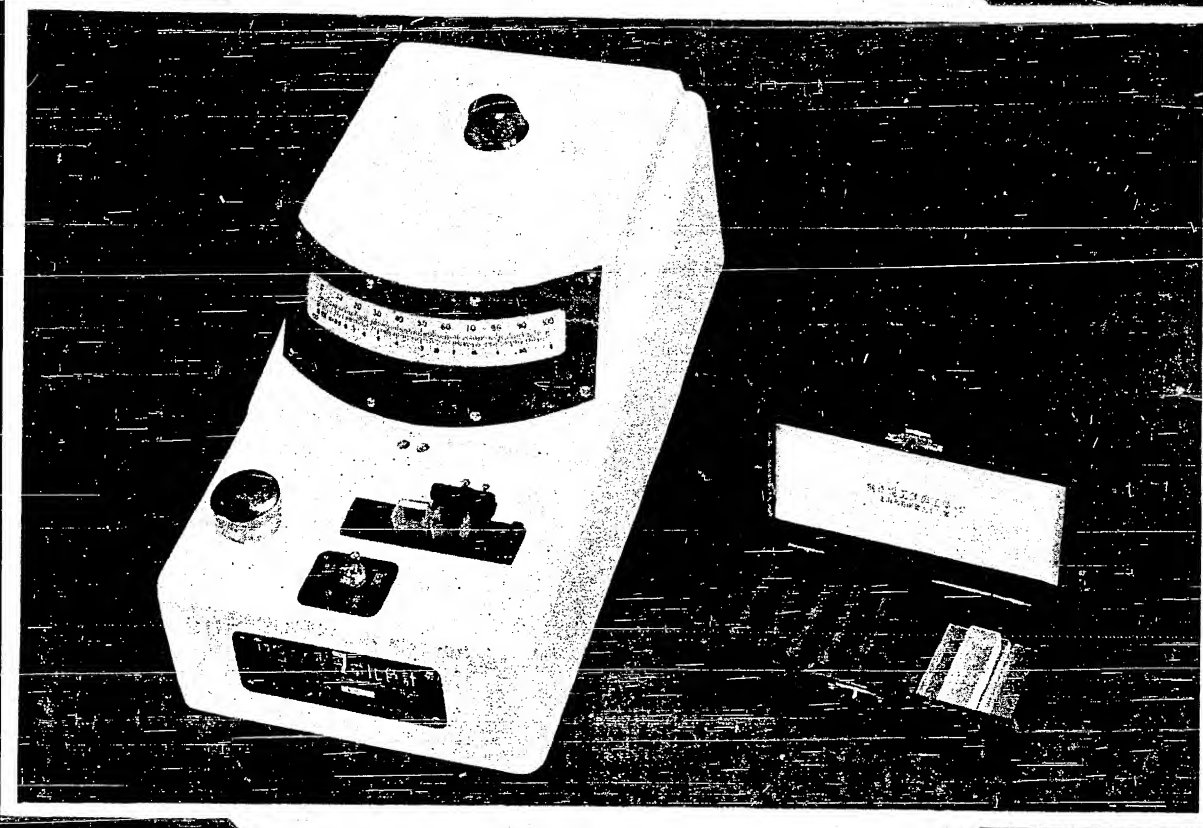CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

PHOTO ELECTRIC COLORIMETER MODEL 93B

The Photo Electric Colorimeter Model 93B embodies the most recent improvements in design for efficiency, compactness and economy. This instrument is extremely simple to use and is easily standardized by reading the transmittancies of a number of solution of predetermined concentrations prepared in accordance with a specific chemical procedure. These readings, when plotted on a semilogarithmic scale against the known concentration on a linear scale, yield an analytical (calibration) curve or standard against which the unknown solutions of that technique may be compared.

This instrument is designed to work either from a 6 volt battery or direct from a 220-250 volt a. c. mains supply through a constant voltage transformer. The light passes through an absorption cell with an optical path of 10 mm., which holds the solution under test. Through the selected filter the transmitted light reaches the photo-cell, the output current of which is indicated on the galvanometer.

The galvanometer is calibrated with a linear scale 0-100 underneath which is a logarithmic scale for calculation of optical density. A set of six filters is supplied.

PHOTOELECTRIC TURBIDIMETER MODEL 71B

STAT


PHOTOELECTRIC TURBIDIMETER MODEL 71B

The Photoelectric Turbidimeter Model 71B is a photoelectric instrument for measuring the light transmission of liquids in Nessler tubes, particularly suited for the measurement of faint turbidity in the field of water testing. It is also suited for the measurement of samples of pale-colour in all colorimetric analysis and colour standardization and grading of light coloured liquid products.

The light source of the instrument can be operated either from the 105-125 volt 60 cycle A.C. power line or from a 6-8 volt storage battery. The light projected vertically upward by a concentrated filament lamp. It is collimated by a lens to form a parallel beam, and passes through the bottom of the tube and the liquid through a colour filter. The light strikes a photocell and is registered by a spotlight galvanometer of high sensitivity. Two scales are provided, the one indicating percent transmission, the other optical density units.

The Reference Standard for turbidity tests, is a suspension of 1000 parts per million of silica (Fullers earth or diatomaceous earth) in distilled water. This suspension of 1000 ppm. has been assigned the value "Turbidity 1000". To calibrate the instrument for turbidity, the 1000 ppm. reference standard is to be diluted with distilled water to furnish turbidities of 10, 20, 40, 60 80 and 100 ppm. These diluted turbidity standards are then read in the instrument, using the blue filter against distilled water as the blank. By plotting the density readings against ppm. a curve will be obtained. Once this calibration curve has been established, in routine test, it is only necessary to measure the sample with blue filter after which the curve furnished the Turbidity value.

Thermistorized Thermometer

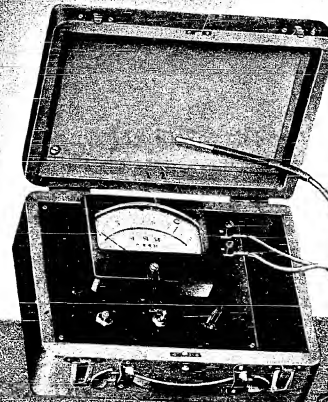CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
SHANGHAI BRANCH
IMPORTERS & EXPORTERS

THERMISTORIZED THERMOMETER

(For Hospital and Medical Use)

The Thermistorized Thermometer can be considered as sensational novelty for hospitals and medical uses. This thermometer is fitted with an electric torch battery (interchangeable) and, therefore, completely independent from the electric circuit.

This thermometer shows immediately the temperature at any part of the human body and gives the following possibilities and advantages to the physician :

1. Immediate indication of the measured temperature, within 1-3 seconds.
2. In order to investigate disturbances of circulation of the human blood, the temperature of the skin can be measured practically at the same time on all parts of the body.
3. Owing to the difference of temperature, it will be possible also to localize inflammations.

PHOTOMETER MODEL 53

中國儀器進口公司
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

MODEL 53 PHOTOMETER

Model 53 Photo Meter is designed specially for measuring intensity of illumination in public auditoriums, schools, hospitals, factories and laboratories. It is a combination of high grade selenium photocell and a precision microammeter. The latter is calibrated to read meter-candle from 0-25000 lux, divided into 3 ranges. A ranging switch is found on the left of the panel. The switch also provides circuit breaking action for safety control to protect the photocell from the strong illumination. The calibration adjusting knob is attached to the right side of the panel. It is easy to adjust the calibration of the meter scale caused by the physical variation of the sensitive microammeter.

The photocell produces a current proportional to the intensity of illumination falling upon it. This instrument may be used for checking either daylight or artificial illumination in laboratories for chemical analysis.

A glass filter sputtered with gold in high vacuum system should cover on the photocell when measuring under strong illumination.

It is simple in operation and gives quick and accurate reading of illumination.

It is compact and light weight with overall dimension of $7\frac{1}{2}$ by $9\frac{1}{2}$ by $4\frac{1}{2}$ inch.

PORTABLE POTENTIOMETER TYPE 703

STAT

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

I APPLICATIONS

Portable potentiometer Type 703 is an instrument for measuring the DC electromotive force. Its range is from 0–71 millivolts.

This portable potentiometer is especially suitable for direct measuring of thermocouples, for comparison against standard thermocouple, or of calibrating temperature indicating millivoltmeters.

II CONSTRUCTION

Circuit diagram of the portable potentiometer is as follows.

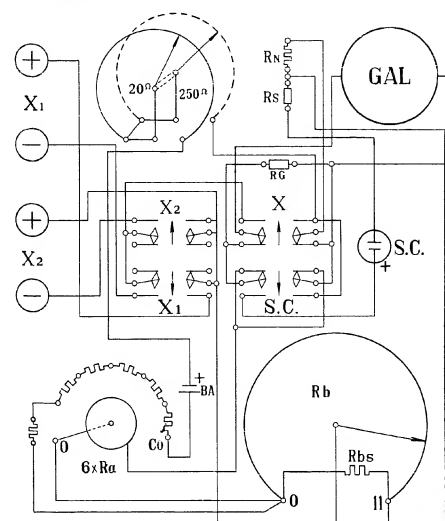
FIG 1

1. The measuring part comprises a range selecting switch and a slide wire. The range switch has six ranges, and the slide wire had a total range of 0-11 MV with 0.1 MV per division.
2. The working current of the potentiometer is 3 MA and is adjusted to this value by comparison with an internal standard cell.
3. The current adjuster comprises a two-ganged rheostat, (fine and coarse rheostats in series). It is so arranged that when turning the knob in one direction the fine and coarse rheostats rotate together and coarse adjustment is obtained; when turning the knob in the opposite direction only the fine rheostat rotates, and thus a finer adjustment is obtained.
4. A pointer galvanometer is used for null indicating.
5. A protecting resistor is inserted in series with the standard cell circuit. When the voltage of the battery falls beyond the adjusting limit of the rheostat, this resistor will protect the S.C. from being damaged by over current.

III TECHNICAL CHARACTERISTICS

- | | |
|----------------------------------------------------------------------------------------------------------|--------------------------------|
| 1. Range of the potentiometer | 0-71 mv |
| 2. Maximum error of the potentiometer at $20^{\circ} \pm 2^{\circ}\text{C}$ and relative humidity 30-80% | ± 0.1 mv |
| 3. Auxiliary error when ambient temperature is $20^{\circ} \pm 10^{\circ}\text{C}$ | less than $\pm 0.2\%$ |
| 4. Specification of the pointer gal. | |
| a/ Internal resistance | 15 ohms |
| b/ Critical resistance | 60 ohms |
| c/ Voltage constant | less than 200 uV/mm |
| 5. Dimension of the case | $192 \times 274 \times 234$ mm |
| 6. Weight of the instrument | 5.5 Kg |

IV COMPLETE EQUIPMENT

The instrument is supplied with following parts:

- | | |
|---------------------------------------------|--------|
| 1. Standard cell (installed in the cabinet) | 1 pc. |
| 2. Dry cells (installed in the cabinet) | 3 pcs. |
| 3. Instruction Manual | 1 pc. |
| 4. Certificate of Test | 1 pc. |

V OPERATION

1. To Measure E.M.F.
 - 1) Release the pointer of the gal by pushing the button backward. Adjust the zero adjuster until the pointer indicates in the zero position.
 - 2) Connect the unknown E.M.F. to the posts X_1 or X_2 . Set the lever switch A to "S.C.". Adjust the current adjuster until the galvanometer indicates zero. The working current is then adjusted to the rated value.
 - 3) Set switch A to "X", switch B to " X_1 " or " X_2 ", balance the potentiometer by turning the range switch and slide wire dial until the gal. indicates zero.
 - 4) Reading on the range switch plus reading on the slide wire dial will be the value measured.
 - 5) Push the lever switch to the middle position when the instrument is not in use. The galvanometer should also be locked.
2. Change battery if the current adjusting rheostat cannot set the gal. to zero position.
3. For calibrating millivoltmeters, connect as shown in Fig. 2

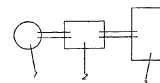
FIG 2

1. Millivoltmeter 2. Variable voltage source 3. Portable potentiometer
4. For convenience of checking ordinary thermocouple against a standard element, connect these two couples to the posts "X₁" and "X₂" respectively. After adjusting the working current, set lever switch A to X, and proceed to measure the couples alternately by setting the switch B to "X₁" and "X₂". Thus any error due to time delay in changing thermocouples may be avoided.
5. The wooden case of the instrument should only be placed in the position supported by the rubber plugs. To place the case in any other position will damage the standard cell inside the instrument.
6. The instrument should be placed in a surrounding temperature of + 5° to + 35°C. relative humidity of 80%. Avoid any direct sunlight.

Note: Mistakes in printing:-

Page 2, line 8:	rheostats for rhostats
" 2, " 14:	circuit for circuit
" 2, " 17:	TECHNICAL for TECHENICAL
" 3, " 6:	Certificate for Cartificate
" 3, " 9:	gal. for gal
" 3, " 22:	use. for use

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMPOT"


BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.L).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

6
MIDGET DUST IMPINGER

MODEL 38B

STAT

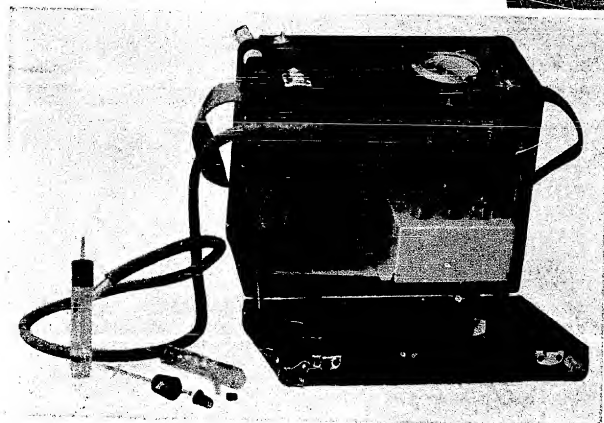中國儀器進口公司
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
IMPORTERS & EXPORTERS

MODEL 38B MIDGET DUST IMPINGER

The model 38B Midget Dust Impinger is a small, compact instrument for standard measurement of atmospheric dust. The importance of the hazard of dust in the industrial plants to health and property is generally recognized. To prevent such disasters an analysis laboratory of dust containment in the air should be well equipped. We are glad to suggest that our Midget Dust Impinger will fulfil all the requirements of dust analysis and brings efficiency to work desired.

This latest development in dust sampling apparatus overcomes the main disadvantages of the regular size impinger—its bulk, weight and power requirements—yet it retains the characteristics which yield reproducible results and make the impinger indispensable for control work and dust surveys.

The pump (Rotary Vacuum Pump) of this impinger, operating in conjunction with a vacuum regulator, maintains an unfluctuating suction regardless of minor variation in the rate or smoothness of cranking.


Isolated locations, confined spaces, the absence of electricity or compressed air are not obstacles to the use of the portable Midget Dust Impinger since it can be carried and operated almost everywhere. This apparatus makes it possible to conduct dust surveys easily and quickly.

The portable Dust Impinger is supplied complete with dust collecting tubes, 9 impinger tubes, rubber, stopper, carrying case and shoulder strap. It is 35.5 cm long, 23 cm wide and 30.5 cm deep, weights 10 kgs.

①

COLONY COUNTER^{STAT}

MODEL 54

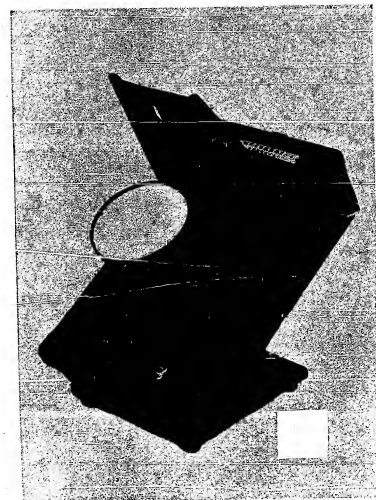中國儀器進口公司
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
SHANGHAI BRANCH
(IMPORTERS & EXPORTERS)

MODEL 38B MIDGET DUST IMPINGER

The model 38B Midget Dust Impinger is a small, compact instrument for standard measurement of atmospheric dust. The importance of the hazard of dust in the industrial plants to health and property is generally recognized. To prevent such disasters an analysis laboratory of dust containment in the air should be well equipped. We are glad to suggest that our Midget Dust Impinger will fulfil all the requirements of dust analysis and brings efficiency to work desired.

This latest development in dust sampling apparatus overcomes the main disadvantages of the regular size impinger—its bulk, weight and power requirements—yet it retains the characteristics which yield reproducible results and make the impinger indispensable for control work and dust surveys.

The pump (Rotary Vacuum Pump) of this impinger, operating in conjunction with a vacuum regulator, maintains an unfluctuating suction regardless of minor variation in the rate or smoothness of cranking.

Isolated locations, confined spaces, the absence of electricity or compressed air are not obstacles to the use of the portable Midget Dust Impinger since it can be carried and operated almost every where. This apparatus makes it possible to conduct dust surveys easily and quickly.

The portable Dust Impinger is supplied complete with dust collecting tubes, 9 impinger tubes, rubber, stopper, carrying case and shoulder strap. It is 35.5 cm long, 23 cm wide and 30.5 cm deep, weights 10 kgs.

⑦

COLONY COUNTER

STAT

MODEL 54

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CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

SHANGHAI BRANCH
(IMPORTERS & EXPORTERS)

MODEL 54 COLONY COUNTER

This Colony counter is an improved instrument and is recommended for Laboratories of Health Department and Hospitals, Water Testing and Milk Testing Laboratories.

It is so designed that this equipment will give maximum testing speed for quick counting of colony in high accuracy. It provides dark field illumination and bright dividing line for making bacteria colony counts.

The case and its stand are made of sheet metal, well built and finely finished on the outside in black wrinkle enamel paint. The angle of the case is adjustable to suit best convenience of the observers by means of adjusting the turning knob of the slot arm.

The 60 watt incandescent bulb is located in the upper end of the Colony Counter and its illumination is intensified by a mirror so placed at the lower end that the light is reflected back to the lower surface of the counting plate, striking it at the same acute angle as the light from the bulb. This principle increases the contrast between the colony and the culture medium. It improves the visibility so that very small colonies can be seen. As a result the observer works efficiently and obtains accurate data quickly.

The $\frac{4}{8}$ inches lens brings the entire dish into view at a magnification of $1\frac{1}{2}$ diameter.

A mechanical counter is attached on the left side of the case, it provides automatic recording by finger pressure. A knurled reset knob is figured to zero.

The equipment is 20 cm. wide, 40 cm. long, 33 cm. deep with two incandescent bulbs for spare.

TYPE 101

STAT

PORTABLE
KELVIN
BRIDGE

INSTRUCTION MANUAL

Portable Kelvin Bridge is a convenient instrument for low-resistance measurement in the range of 0.0001–10 ohms, or for measuring the conductivity of wires.

The circuit diagram of the bridge is as follows:

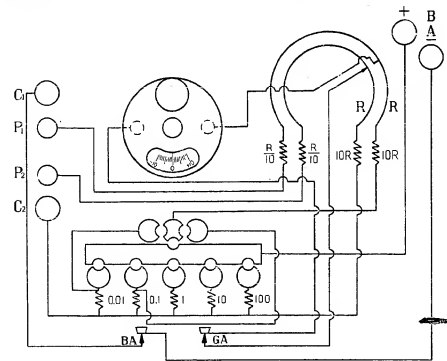
FIG 1

CONSTRUCTION

A knob at right of the panel is used to vary the position of the galvanometer circuit in the double slide wire, and therefore the reading of the dial. The galvanometer is installed at the left of the panel, with a ratio block in the lower part of the meter.

(1) The galvanometer

A zero adjusting cap (7 in fig. 2) is mounted on the upper part of the galvanometer. Care should be taken not to turn violently to avoid damage of the adjuster. A button (9 in fig. 2) is provided to lock the pointer of the galvanometer. To lock the pointer, push the button forward, to unlock, push backward.

(2) Ratio Block

Two plugs are used to change the ratio of the bridge (10 and 11 in fig. 2). Plugs should be pressed firmly to ensure minimum contact resistance.

(3) Slide Wires

The reading dial comprises two slide wires. Turning the knob of the dial will vary the position of the galvanometer in the slide wires.

(4) Keys

Key "G" (6 in fig. 2) may be locked by pressing down and turn clockwise.
Key "B" (5 in fig. 2) cannot be locked.

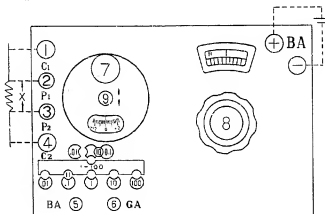
FIG 2

- (5) Weight of the instrument approx 3.4 kg.
Dimension shown as follows:

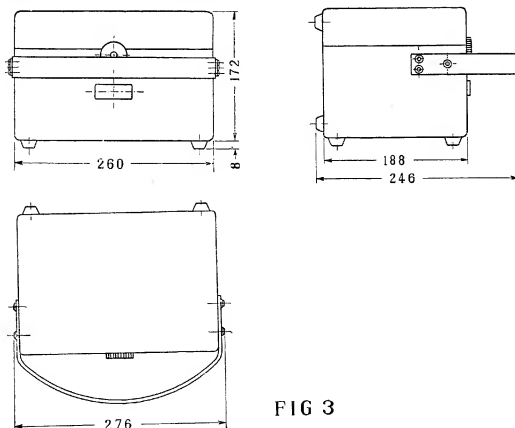
FIG 3

OPERATION

- (1) Resistance lower than one ohm must have four terminals, connect as shown in fig. 2. Posts No. 1 and No. 4 are current terminals, and posts No. 2 and No. 3 are potential terminals. For resistance higher than one ohm, posts No. 1 No. 2, and No. 3 No. 4 may be tied up with a short bus wire respectively, and resistance to be measured is connected to posts No. 2 and No. 3.
- (2) Connect battery to posts "B+" and "B-". Use three 1.5v dry cell (connect in parallel), or one 2v 60Ah storage battery. When measuring low resistance, better use 4v storage battery. Key "B" (5) should not be pressed too long to avoid excessive power consumption.
- (3) Press down key 6 to close the galvanometer circuit. To lock the key, press down and turn clockwise.
- (4) Release the galvanometer lock 9, and turn the zero adjuster 7 until the pointer indicates in the zero position. Set the plugs in the corresponding holes as follows:

Measuring Range	Position of Plug 10	Position of Plug 11
1—10 Ω	1—100	100
0.1—1 Ω	1—100	10
0.01—0.1 Ω	1—100	1
0.001—0.01 Ω	0.1	0.1
0.0001—0.001 Ω	0.01	0.01

- (5) Press down key 5, and balance the bridge by turning the dial 8. The dial reading multiplied by the value indicated at plug 11 will be the result of the unknown resistance.

EXAMPLE

The dial reading 0.025

The plug in the 0.1 position

$\therefore R_X = 0.025 \times 0.1 = 0.0025 \text{ ohm}$

- (6) The null indicator of the bridge is a sensitive suspension type galvanometer avoid excessive vibration. Push forward the button 9 to lock the pointer when not in use.

- (7) To measure the conductivity of copper wire, use following formula:

$$\text{Conductivity} = R_t / R_X \times 100\%$$

R_X = Resistance of copper wire to be measured at $t^\circ\text{C}$

R_t = Resistance of standard copper wire of same diameter at $t^\circ\text{C}$

R_t may be calculated as follows:

$$R_t = R_{20} [1 + 0.003(t - 20)]$$

where R_{20} is the resistance of standard copper wire of same diameter at 20°C .

- (8) The instrument should be kept in dry place, and should avoid any harmful atmosphere.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMPOR"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.I).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

PORTABLE
WHEATSTONE
BRIDGE

STAT

TYPE 102

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CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
IMPORTERS & EXPORTERS

I APPLICATIONS

Portable wheatstone Bridge is designed for general measurements of resistance from 1 to 107 ohms, and for locating faults in electrical cables.

The circuit diagram is as follows:

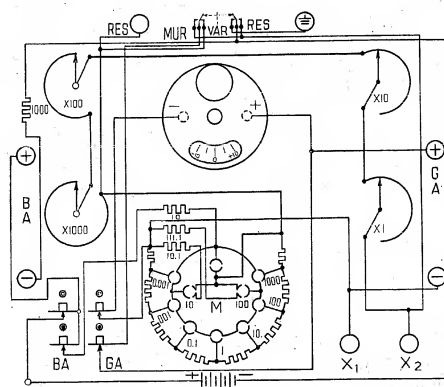
FIG 1

II CONSTRUCTION

The panel of the instrument is made from high quality bakelite, therefore good insulation is ensured. Four knobs for varying the resistance units are mounted on the panel, and a lever switch is placed at the upper part of the panel. The lever switch has three positions i.e. RES, MUR, VAR. A pointer galvanometer is situated in the middle of the panel. Below it is the ratio arm plug. The bridge has eight binding posts, and four keys are provided for gal. and battery.

- 1) Galvanometer. A zero adjusting cap is mounted on the galvanometer. Turning the cap will adjust the zero position of the pointer. Do not turn violently to avoid damage to the pointer. To lock the pointer, when not in use, push forward the button.
- 2) Switch Dial. Four bakelite knobs are mounted on the panel to vary the resistance of the measuring arm.
- 3) Ratio Arm. The ratio arm has 10 plug holes. Three holes for Murray loop test are surrounded by 7 holes for resistance measurements.
- 4) Lever Switch. The switch is mounted on the upper part of the panel. It has three positions, i.e. RES, MUR, VAR.
- 5) Keys: Four keys are situated at the lower left corner of the panel. Keys 3 and 7 (Fig. 2) are for gal. and 2 and 6 are for battery.


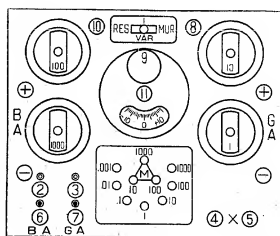
FIG 2

- 6) Case: A small compartment for installing dry cells is provided at the lower part of the case. A thick aluminum plate is used as a cover for the compartment.
- 7) Weight of the instrument: 3.9 Kg. Dimensions: as in Fig. 3.

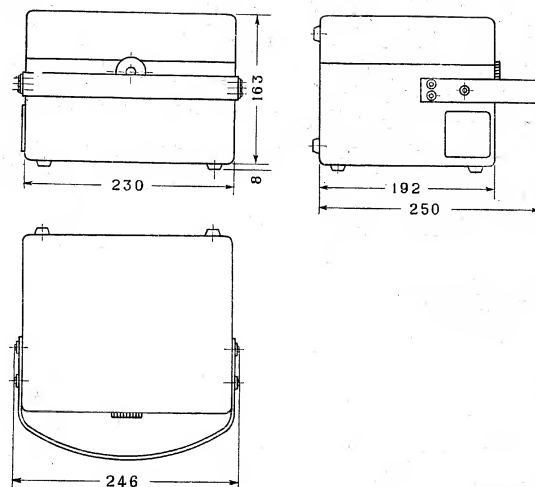III OPERATION

I Resistance Measurements

A. Resistance ranging from 1—10⁵ ohms.

- (1) Push the lever switch 1 towards "RES". Release the gal. lock by pushing the button downward. Adjust the pointer to zero position.
- (2) Connect the unknown resistance R_x to posts "X₁" & "X₂".
- (3) Lock the keys 2 & 3 (Fig. 2)
- (4) Set the plug to the corresponding ratio as shown in the table below

Unknown Resistance	Ratio
<10Ω	0.001
10—100	0.01

100--1000	0.1
1000--10000	1
10000--100000	10

Table 1

Internal Gal. & Bat.	External Gal. Internal Bat.	External Bat. Internal Gal.	External Gal & Bat
Press Gal	Press Gal	Press Gal	Press Gal

(5) Balance the bridge by varying the rheostat and at the same time press down keys 6 & 7. Positions of the keys for connecting external battery and galvanometer are shown in table 1.

(6) Let the value of the measuring arm is R, and ratio A, the unknown resistance R_X is

$$R_X = A \times R$$

B. Resistance ranging from 10^5 — 10^7 ohms.

- (1) Connect as shown in Fig. 4.
- (2) For keys and gal. same as in A (1)
- (3) For external battery connection

release key 2 and close k 1. for external gal., release key 3 and lock the pointer of the gal. in the bridge, then close key 2.

- (4) The ratio setting, voltage applied, and sensitivity of the gal. should be as follows:

Unknown Res.	Ratio	External Voltage	Sensitivity of Gal.
10^5 — 10^6 ohms	100	15	100mm/uA
10^6 — 10^7 ohms	1000	25	> 100mm/uA

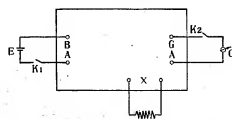
FIG 4

II Used as a resistance box

Release all the keys, and posts 5 & 10 (Fig 2) will be terminals of an individual four decade resistance box.

III Locating fault cables

Fault cables such as breakdown of insulation, breaking of wire, may be located by this instrument.

A. Murray Loop Test

Murray loop test is used for test of ground of a single wire cable. Connect as shown in Fig. 5.

- (1) Join the faulty wire to a good one at the distant end of the cable; connect faulty wire to X_2 and good wire to X_1 . Connect post 8 to ground. Set ratio plug to M1000. Push the lever switch at MUR-Vary rheostat to balance the bridge, if a satisfactory balance cannot be obtained with M1000 change to M100 or M10.

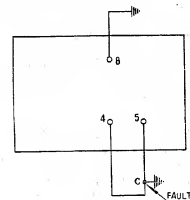
FIG 5

- (2) Measure the total resistance of the wire and call it L.
- (3) If the total length of the cable is L' , then Resistance R_X between 5 and C may be calculated from the following:

$$R_X = \frac{L \cdot R}{M + R} \text{ ohms}$$

$$\text{Length of X} = \frac{2L \cdot R}{M + R} \text{ ohms}$$

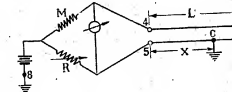
FIG 6

B. Varley Loop Test

- (1) Join the faulty and good wires at the distant end of the cable; connect the faulty wire to X_2 and the good wire to X_1 . Connect post 8 to ground. Set lever switch at VAR. Set ratio plug A at 0.1 or 0.01 (Fig 7) and vary rheostat R to balance the bridge.

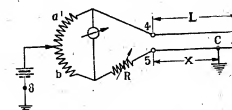
FIG 7

- (2) Measure the total resistance of the wire and call it L .
 (3) If the total length of the cable is L , then resistance R_x between 5 and C may be calculated from the following formula:

$$R_x = \frac{L - AR}{A + 1}$$

If the length of the cable is L , then length of X will be:

$$X = 2L' \frac{R_x}{L}$$

C. Three Varley Method

This method is used to locate the fault of three wire cable.

- (1) Join the three wires at distant end of the cable; Connect faulty wire to X_2 , one good wire to X_1 , the other good wire to post 8. Push lever switch to VAR. Set the ratio plug to an appropriate value A . Balance the bridge and obtain reading R_1 (Fig 8)

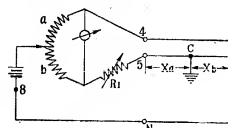
FIG 8

- (2) Connect post 8 to ground instead of good wire. Balance the bridge with the original ratio, and obtain reading R_2 (Fig 9) Resistance of X_b is

$$X_b = \frac{A}{A+1} (R_2 - R_1)$$

- (3) Push lever switch to "RES". Balance the bridge with the original ratio, and obtain reading R_3 (Fig 10)

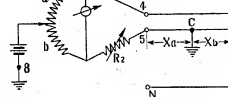
FIG 9

Resistance of X_a is

$$X_a = \frac{A}{A+1} (R_3 - R_2)$$

- (4) Let the length of the cable be L' . distance from 5 to C may be calculated as follows:

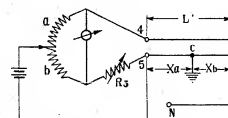
FIG 10

$$X = L' \frac{X_a}{X_a + X_b}$$

D. Open Location-Quaded Cable

This method requires external low frequency source (<100 c/s) and headphone for open location of quaded cable.

- (1) Connect the open wire to post X_2 and a good wire of equal gauge and length in the same cable to X_1 , and strap them together at the distant end of the cable. Join the mates of these wires at both end. All wires must be free from grounds and crosses. Connect one terminal of source of A.C. to the -BA post and the other terminal to the near end of the joined mates. Set ratio plug at M1000, M100, or M10. Set lever switch at MUR. Depress keys-6 & 7 and balance bridge by varying the rheostat for silence or minimum tone in headphone. (Fig. 12)

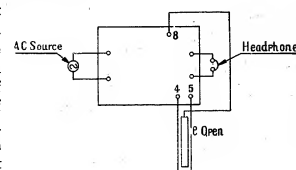
FIG 11

If the length of the cable is L' , then distance from 5 to open location C may be calculated as follows:

$$X = \frac{2L'M}{M+R}$$

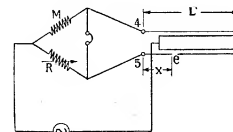
FIG 12

E. Open Location-Pairs

- (1) Connect low frequency source (<100 c/s) to the posts for external source. Set the lever switch to MUR. Connect the open wire to X_1 and ground its distant end. Ground the near end of its mate with distant end open. Connect X_2 to one terminal of a condenser with the other terminal grounded. Arrange bridge as in preceding description, but ground outer terminal of source, Balance the bridge and obtain reading R_1 (Fig 13)

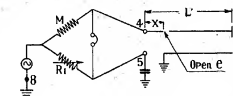
FIG 13

- (2) Disconnect near end of faulty wire from X_1 and ground it. Disconnect

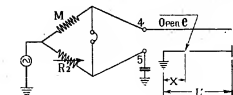
FIG 14

near end of good wire from ground and connect it to X_1 . Balance the bridge and note reading as R_2 (Fig. 14)

- (3) If the length of the cable is L , then length X may be calculated as follows:

$$X = \frac{L R_1}{R_2}$$

IV BATTERY MAINTAINANCE

To install battery, remove the aluminum plate at the left side of the case, and insert the cells. Change battery at least once every six months. (Fig15)

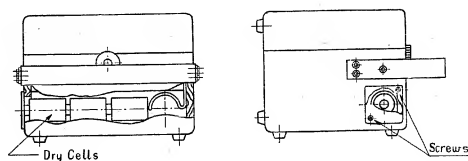
FIG 15

V CAUTION

- (1) For measuring resistance containing inductance, key 7 should be released. Depress key 6 and then key 7.
- (2) Release key 3 when using external gal. The pointer of the gal. in the instrument should be locked.
- (3) Release key 2 for external battery.
- (4) Lock the pointer of galvanometer when not in use.
- (5) The instrument should be placed horizontally when measuring.
- (6) Battery should be changed periodically to avoid damage to the instrument.
- (7) Excessive vibration should be avoided.
- (8) The instrument should be kept in a cool dry place. Keep the panel to ensure good insulation.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou. Hsi Chiao, Peking, China.
Cable Address: "INSTRIMPORT"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E. 1).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 2 West Bund
Cable Address: "INSTRIMP"


A G E N T:

HONGKONG: China Resources Co.
Bank of China Building 11/12 floor, Hongkong
Cable Address: "CIRECO HONGKONG"

**PRECISION
WHEATSTONE
BRIDGE**

"STAT

TYPE 103

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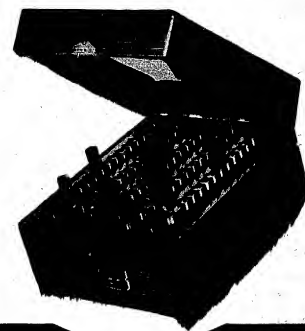CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

RECEIVED 8 12 1964

I. APPLICATION

Type 103 Precision Wheatstone Bridge is intended for direct current resistance measurements of the highest precision in laboratory. The measuring range is from 10^{-5} to 10^8 ohms, readable to five digits.

Two binding posts are installed on the panel from the measuring arm so that it can be used as an individual five decade plug type resistance box.

II. CONSTRUCTION


Type 103 Precision Wheatstone Bridge consists of two parts, the measuring arm and the ratio arms. (Fig. 1)

1. The Measuring Arm

The measuring arm has five decades of ten resistors each, respectively, 0.1, 1, 10, 100, and 1,000 ohms; limit of error for 0.1-ohm resistors $\pm 0.1\%$, for all others $\pm 0.05\%$. The bridge blocks have five plugs which can be pressed firmly to ensure small contact resistance.

2. The Ratio Arm

Ratio arms have two each of 1, 10, 100, 1,000, 10,000-ohms resistors; limit of error in ratio resistors $\pm 0.02\%$. Changing of the ratio is obtained simply by changing the positions of the two plugs. Plugs for exchanging the positions of


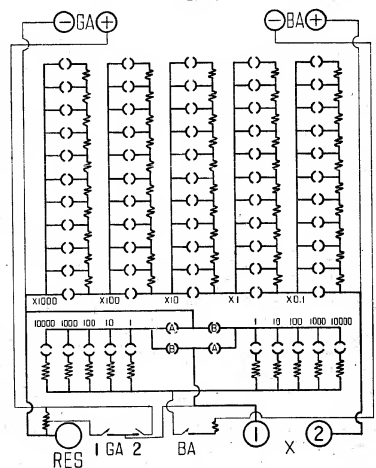
CIRCUIT DIAGRAM

of the ratio arms are also provided. For precise measurement, first measure one value with arms set at "A", then exchanging the position of the arms and measure again, the mean value will eliminate any possible error due to the contact and wiring resistance of the ratio arms.

3. A shunt is provided in the galvanometer circuit. When measuring, first press down key GA—1 and balance the bridge, then release GA—1, and press down key GA—2, the balanced value will be the desired value.
4. A 10—ohm resistor is installed in series with the battery circuit to protect the ratio and measuring arms from being damaged by over current.

III GENERAL CHARACTERISTICS

1. Ratio arms 2 (1+10+100+1000+10000) ohms, max. error $\pm 0.02\%$
2. Measuring arm 10 (0.1+1+10+100+1000) ohms, max. error $\pm 0.05\%$ for one ohm and higher, $\pm 0.1\%$ for lower than 1—ohm.
3. Measuring range 10^{-5} — 10^8 ohms.

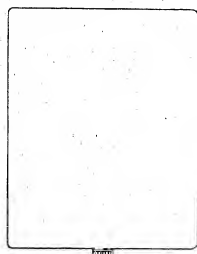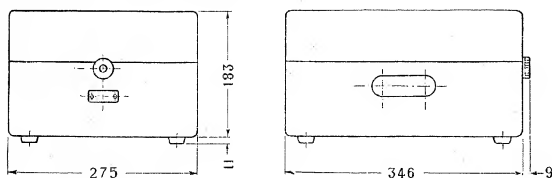
EXTERNAL VIEW

4. Maximum resistance of the measuring arm, 11,111 ohms.
5. Weight 9.8 kg
6. Dimension 355 x 275 x 194 mm (Fig. 2)

IV OPERATION

A. Resistance measurement

1. Connect as in Fig. 3

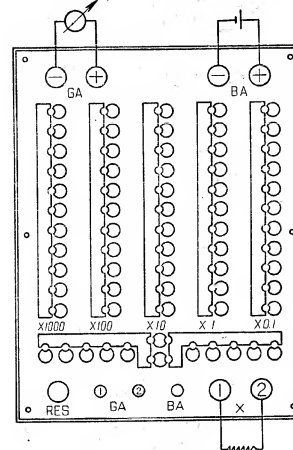CONNECTING DIAGRAM

2. Estimate the unknown resistance and choose the corresponding ratio:

Unknown Res.	Ratio	Voltage	Unknown Res.	Ratio	Voltage
1	1/10000		10^3 — 10^4	1	8
1—10	1/1000	2	10^4 — 10^5	10	10
10—100	1/100	4	10^5 — 10^6	100	20
100—1000	1/10	6	10^6 — 10^7	1000	
			10^7 — 10^8	10000	

3. Set the plugs in the measuring arm to some appropriate value so as to make the unknown resistance equal to ratio X measuring value.
4. Balance the bridge by changing the plug position until the galvanometer does not deflect when GA key is depressed. Balance first with GA-1, then with GA-2. The unknown resistance equals to the measuring plug-setting multiplied by the ratio arm plug-setting.
5. For precise measurement, it is desirable to exchange the ratio arms in order to have a mean reading. For example, with arm choosing plugs in position "A" and ratio plugs right in 1000 and left in 100, we have a certain reading, then change the arm choosing plugs to position "B", and ratio plugs right in 100 and left in 1000, add the reading to the first reading and divide by two, we obtain the mean value.


B. Used as a plug type resistance box.

1. Connect to the posts "RES" and "X2" as shown in the diagram.
2. Set the plugs in the desired position, for example 1643 ohms, we have the plugs in hole No. 1 of the thousandth decade, No. 6 of the hundredth decade, No. 4 of the tenth decade and No. 3 of the unit decade. The plug in the final decade being put in 0 position. The resulting resistance will be 1643.0 ohms.

Note: The maximum rated power for each coil is $\frac{1}{4}$ watt.

C. Auxiliary Equipment required and should be purchased separately.

1. Power source of 2-20 volts (Battery or rectifier)
2. Reflecting galvanometer with sensitivity not lower than 100 mm/microamp.

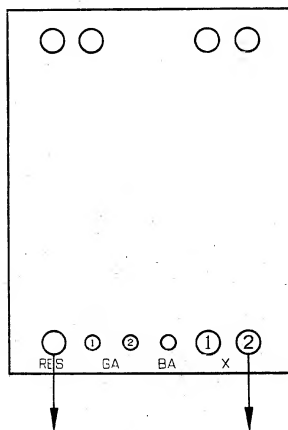V PACKING

The cartoon contains the following apparatus and accessories:

- | | |
|-----------------------------------------|---|
| 1. Type 103 Precision Wheatstone Bridge | 1 |
| 2. Brass Plugs | 9 |
| 3. Instruction Manual | 1 |
| 4. Certificate | 1 |

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMPORT"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.1).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

**A.C. SINGLE PHASE^{STAT}
WATTHOUR METER**

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**CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)**

A. C. SINGLE PHASE WATTHOUR METER

The S-1 model watthour meter is a single two-wire meter of our latest design. With skilled experience and improvement the S-1 model is constructed in a beautiful totally enclosed case, so as to prevent moisture and dust from entering into the inner parts and to assure a long service.

1. USES:— The S-1 model is a single phase two-wire induction type for measuring A.C. current 50 cycles and gives accurate reading of KILOWATT HOUR in accordance with consumption and period of time in operation.

2. CONSTRUCTION:— S-1 watthour meter weighs about 1½ kg. Back plate and front casing are made of sheet steel with a piece of glass fitted in the cover to make the movement of the revolving disc visible. All the parts and accessories are mounted on the inside frame. Electro-magnetic parts, including potential coil, current coil, temperature compensating coil, overload compensator, 90° lagging coil, 90° adjuster and light load adjuster, are made from silicon steel laminations and are fitted on the left side of the frame. The 90° adjuster is of resistance wire connected to the temperature compensation coil for adjusting the phase position of the magnetic flux produced by the current coil. The light load adjuster is a screw rod, which can be turned at both left and right ends for adjusting the speed of the revolving disc at light loads. The revolving parts are fitted on the axle at the front in the middle with the disc extended to the gap of the electro-magnet and permanent magnet. The disc will be rotated by current induced from the current in potential coil and current coil. The permanent magnet is fitted on the mounting frame at the middle of the right part. Its magnetic flux, passing through the disc, creates blocking effect enabling the speed of the revolving disc in direct proportion to the current and voltage of the load. By the side of the permanent magnet, there is a full load adjusting screw to shift the position of the permanent magnet, so to vary the magnetic force, and the speed of the revolving disc at full load can be adjusted accordingly. These features render adjustment easier than that

for meters of other types. On the top of the mounting frame, there is the registering mechanism, the gear wheel of which meshes that of the wheel on the disc axle, the readings of kilowatt-hour consumed can thus be obtained from the cyclo dials with 5 figures. Terminal box is below the casing with four openings. Upper bearing of the axle is of needle sleeve type, while the bottom bearing of the axle is of steel ball and jewel with spring to reduce the wear and increase the sensitivity. The revolving disc has a diameter of 90 mm and weighs 18.8 grams.


3. RATINGS:—

VOLTAGE AND CURRENT RATINGS

Rated Voltage	Rated Current (Amp.)												
110 v.	1½	2½	3	5	10	15	20	25	30	40	50	75	100
220 v.	1½	2½	3	5	10	15	20	25	30	40	50	75	100

Current transformer will be supplied for current over 100 Amp.

4. DIAGRAM OF CONNECTION:—

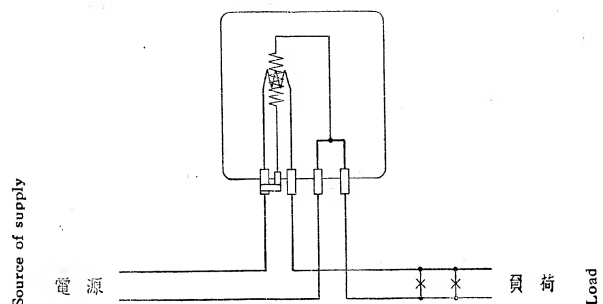5. CHARACTERISTICS:—

1. ACCURACY:

Table (I)

Rated current	Power factor	Permissible error
10 - 125%	1	$\pm 2\%$
5%	1	$\pm 3\%$
10 - 125%	0.5 (lagging)	$\pm 2\%$
10%	0.5 (lagging)	$\pm 3.5\%$

2.

Table (II)

Rated current	Power factor	Permissible error
10 - 125%	1	$\pm 2\%$
10 - 125%	0.5 (lagging)	$\pm 3\%$

3. STARTING CURRENT:— 0.5% — 1% of rated current.

4. CRAWLING:— There is no distinct movement even though the voltage is increased by 110%, and if there is any, it will not exceed one revolution.

5. EFFECT OF VOLTAGE:— The permissible error will not exceed 1% when the applied voltage is 10% over or under the rated voltage.

6. EFFECT OF TEMPERATURE:— When temperature is varied between 0°C and 40°C, the permissible error will not exceed $\pm 2\%$ at unity power factor and $\pm 3\%$ at power factor of 0.5 lagging.

7. INSULATION:— The insulating resistance is 100 megohms between coil and steel case.

8. POTENTIAL TEST:— The meter can stand a potential test at 2000 volts A.C. 50 cycles for a duration of one minute.

9. POWER LOSS:— 2.4 VA for voltage coil.

10. FULL LOAD SPEED:— 22 revolutions per minute.

6. ADJUSTMENT:— Every meter is properly and carefully tested and calibrated before leaving factory. If re-calibration is required due to long use and damage the following adjustment can be made in accordance with the permissible error mentioned above.

1. FULL LOAD ADJUSTMENT:— Connect a load of 100% rated current at unity power factor, and adjust the position of the permanent magnet with the screw at its side by checking with a standard meter.

2. 90° ADJUSTMENT:— Connect a load of 100% rated current at 0.5 power factor and adjust the length of the resistance coil by checking with a standard meter.

3. LIGHT LOAD ADJUSTMENT:— Connect a load of 10% rated current at unity power factor and adjust the revolving rod either to left or right by checking with a standard meter.

7. INSTALLATION:— If the cover of the meter has been opened before installation, it is necessary to clear off all dust from the interior and then close tightly and seal the meter. The meter should be installed at a place where ventilation, temperature and humidity is most suitable.

METHOD OF INSTALLATION:

1. After hanging the meter on the meter board, the cover of the terminal box is removed. The meter is secured to the board by two wood screws one for each lug hole at two sides of the terminal box. Care must be taken that the axis of the revolving disc be perpendicular to the ground level, otherwise accuracy will be effected owing to tilt of the disc.



2. Power lines and load lines are to be connected in accordance with the diagram of connection. Observe whether the revolution of the disc is normal (left to right direction).

3. If the motion is in order, it is necessary to see whether there is any crawl after the load is disconnected.

4. If the crawling movement is within the limit allowed, then the load lines can be connected and have the cover of terminal box closed and sealed.

5. It is important that the load should be in accordance with the rated current of the meter. Any overload will cause a quick run of the revolving disc resulting in the rise of temperature in current coil which will endanger the meter or shorten its life.

8. OVERALL DIMENSIONS:—

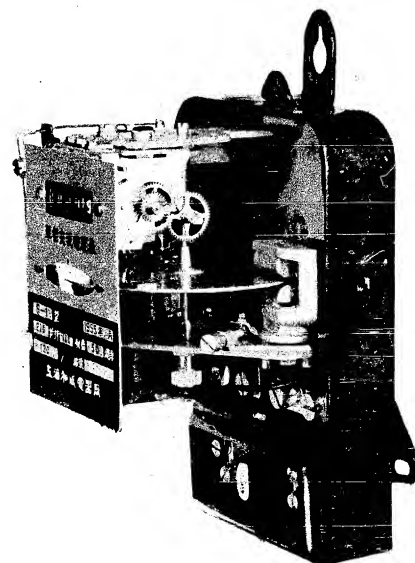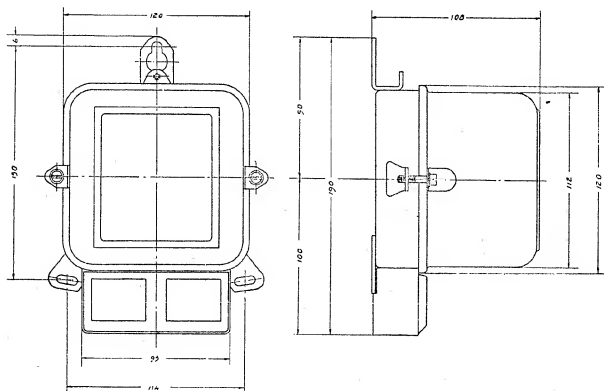CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMP"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.1).
Cable Address: "INSTRIMP"

TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"

CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

STAT

"Chungtien"
HEADPHONE

中國儀器進口公司
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
IMPORTERS & EXPORTERS

CHUNGTIEN HEADPHONES

I MODELS:

2000 Ohm set. 3000 Ohm set. 4000 Ohm set.

II TECHNICAL DATA:

1. The D.C. resistance of the double headphone and turns of coils are as follows:

D. C. Resistance	Resistance of each coil	Turns of coil (Winding)	Wire Gauge
2000	500	3600 turns	SWG No. 46
3000	750	3780 turns	SWG No. 47
4000	1000	5500 turns	SWG No. 47

Note: The above figures are allowed a tolerance limit of 10%.

2. The insulation resistance between the coils and cores and between the cores and bakelite receiver pieces is over 20 Megohms.

3. Permanent magnetic alloy is made of ALNI III magnet alloy.

4. Test of sensitivity. Tone quality is clear and distinct; as a receiver for the crystal radio, its sensitivity is higher than average local product. For simply testing, put one terminal of cord in touch with tongue and let the other terminal of the cord contact slightly with the wrist watch, a clear click will be heard from the receiver.

5. The iron cores are made of magnet iron and properly annealed.

6. The diaphragm is made of 0.17 mm. iron plate also properly annealed.

7. The enamelled copper wire is alloyed less than 10 pin holes per meter.

III This double headphone (or crystal radio headreceiver) is generally accepted by the users for its high sensitivity and easy handling. For maintaining quality, the following vital parts have been specially treated in the course of manufacturing.

1. The diaphragm of the receiver attained high sensitivity by annealing.

2. The raw material for cores is special magnet iron and the annealing process renders the magnetic residual to a minimum.

3. No casting gate is allowed for reproducing the magnetic alloy. The magnetic power is excellent. The rate of degaussing is negligible.

4. The gap between the core and the diaphragm is standardized by precisely processing; the hearing quality is guaranteed.

IV PRECAUTIONS:

1. To prevent the thin copper wire of the coil from oxidation, always keep the phone in a dry place.

2. Do not remove the diaphragms constantly, it effects degaussing.

3. Keep the diaphragm in perfect flatness. Any bending of it brings about defect in sounding quality.

4. Do not connect this phone with a D.C. circuit. A connection regardless of the direction of the current may reduce the magnetic power.

V USES:

This headphone is designed for the crystal radio receiver which is widely adopted in rural areas where electricity is not available. It is also used as a substitution for wire speaker in hospitals and sanitariums where quietness is desired.

This headphone is also available for measuring instruments, wire communications, telegraph coding and other purposes.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMP"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.I.).
Cable Address: "INSTRIMP"

TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"


CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

AGENT: —

HONGKONG: CHINA RESOURCES CO.,
12 Floor, Bank of China Bldg.,
De Voux Road Central,

STAT

“Culture”

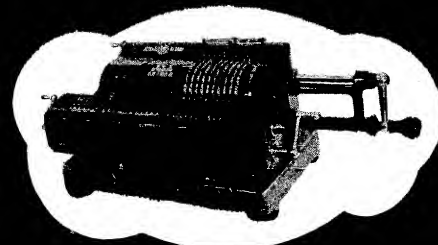CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

"CULTURE" Calculating Machine

Widely employed and well reputed in Commercial, Industrial, Educational and Scientific Circles.

Capacity: $10 \times 10 \times 20$

1. Dependable and Easy operating.
2. Machine parts heat-treated withstanding friction, Chrome-plated and rust proof.
3. Speed: 300 cycles per minute.

Overall:		Packing:
Length:	36 cm.	49 cm.
Width:	19 cm.	30.5 cm.
Height:	13.5 cm.	30 cm.
Nett Wt.	8 kgs.	Gr. Wt. 17 kgs.

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CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION


HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMP"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.1).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

OUR HAND OPERATED
CALCULATING MACHINES
WILL GIVE YOU ALL ANSWERS

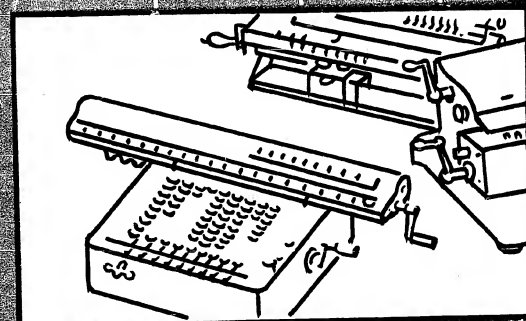TWIN WHEEL

REGISTERED

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

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

(IMPORTERS & EXPORTERS)

"FLYING HORSE"


Calculating Machine Capacity: $10 \times 10 \times 20$

Our "Fly-Horse" calculator is of high efficiency and durable dependability but offered at very low price. Each machine, while reflecting our expert engineering and skill workmanship, is subjected to rigid test and inspection before leaving factory, therefore ensures complete satisfaction.

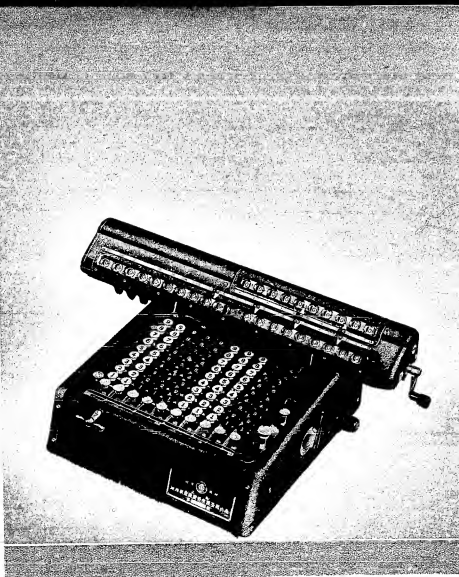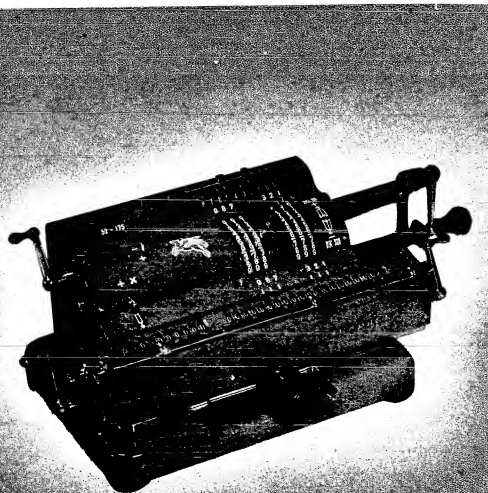"TUNG YUNG" Calculating Machine Capacity: $10 \times 10 \times 20$

"Tung Yung" calculating machine is one of the best high speed hand model calculating machine with full keyboard, easy use, four fundamental calculations, namely: Addition, Subtraction, Multiplication, Division. Besides, it is quite fit for making rapidly several special calculations like: extraction of square roots, elevation to powers of various exponents, cubatures, discounts, calculations with complex numeral or metrical systems etc.

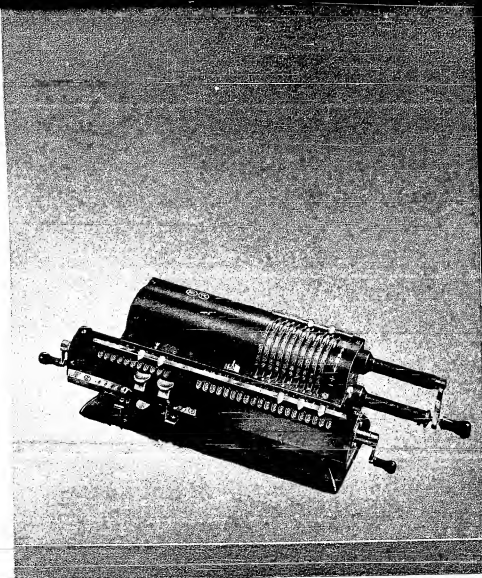"TWIN WHEEL" Calculating Machine

Capacity:
Large register 20 dials
Small register 11 dials
Setting board 10 levers


Adding, multiplying, subtracting and dividing crank type calculating machine. Fast, smooth and easy operation. Handy and light but sturdy efficient machine. Characteristic for this model is its great capacity which is sufficient even for the largest calculating operations that may have to be dealt with.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
 HEAD OFFICE:
 Er-Li-Kou, Hsi Chiao, Peking, China.
 Cable Address: "INSTRIMP"
 BRANCH OFFICES:
 SHANGHAI: 27 Chungshan Road (E.1)
 Cable Address: "INSTRIMP"
 TIENTSIN: 171 Kien Shieh Road.
 Cable Address: "INSTRIMP"
 CANTON: 25 Tai Ping Road, S.
 Cable Address: "INSTRIMP"

STAT

*Cotton Drawing Apparatus.
 Cotton Sorter.*

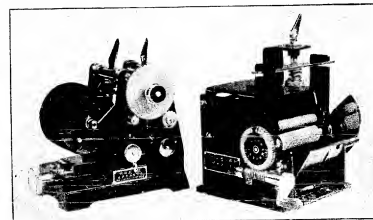中國儀器進口公司
 CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
 (IMPORTERS & EXPORTERS)

HOW TO USE THE DRAWING APPARATUS:

1. Divide the Small Average Test Specimen Cotton into four equal parts.
2. Turn Screw (17) to adjust the distance or space between the two pairs of Rollers to obtain a space of about 3-5 m.m. in excess of the Length Standard under test.
3. Set on Levers (7) to tighten the two pairs of Rollers.
4. Feed the said four equal parts one after another separately into Cotton-feeding Rollers (1) (2); turn Handle (8) counter-clockwise.
5. Take away the foreign matter and neps from the Specimen Cotton on Velvet Roller (12).
6. Take off Specimen from Velvet Roller (12); pass it through the Drawing Apparatus for three to five times (not less than three times for cotton grade 3 and up—not less than five times for cotton grade 4 and down) and make it into four Slivers.
7. Part each of the four Slivers into two equal lengths, discarding one half of each and keeping the other four halves. Combine each two of the remaining four halves into one and pass them again through the Drawing Apparatus for a same number of times as aforesaid to produce two Slivers.

8. Cut each of the two Slivers into two equal halves, again discarding one half of each and keeping the other two halves. Combine the remaining two halves into one and pass it again through the Drawing Apparatus for a same number of times as aforesaid, thus producing the Test Specimen Sliver.
9. Select out, lengthwise, 0.1 to 0.12 gram from the Test Specimen Sliver. Manipulate the fibres carefully by fingers and cast off the dust and other impurities. Again, pass it through the Drawing Apparatus to produce the Ultimate Cotton Sliver.

Cotton Drawing Apparatus


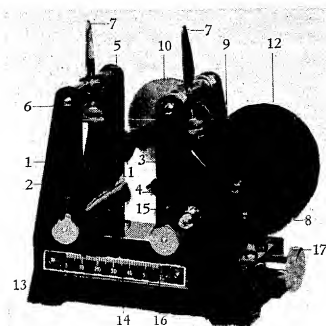
Diagram showing the Process of making Test Specimen Sliver, Ultimate Sliver and Cotton Tuft.

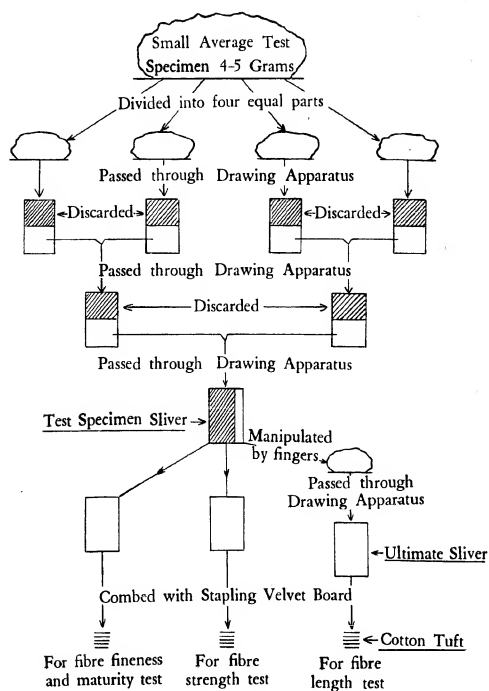HOW TO USE THE COTTON SORTER:

Open the Lid (9) of Cotton Sorter. Set the pointer on the Revolvable Dial (4) at degree 9. Turn downward the wooden part hinged upon Stapling Velvet Board and, with the aid of Fibre Clip No. 1, carefully clip up the cotton tuft which has already been properly made and now placed on Stapling Velvet Board, then place it on the lower roller of Cotton Sorter. When placing the cotton tuft on the roller, the lower adjustable tip of Fibre Clip No. 1 should firmly and properly touch the edge of Sliding Plate (14), and the Clip should also rest upon Support (20) in such a position as to keep the cotton tuft in a horizontal position.


Now, put on the Lid (9) and tighten the Torsion Spring (8). Simultaneously, check the pressure of the Spring with a proper instrument to obtain a pressure of 7000 grams. After taking off the Clip, the even end of the cotton tuft should be closely in line with Sliding Plate (14); thus the exposed part of cotton tuft will be 9 m.m. in length. Now, turn down Sliding Plate (14) in order not to interfere with the work of clipping the short fibres. Then, turn Handle (6) counter-clockwise one round so that the Roller will deliver 1 m.m. of tuft. At this time, the part of the fibres free from the grip of Rollers, within 10 m.m. in length, will be the First Batch.

Draw out this First Batch of fibres with Fibre Clip No. 2 and place it on the Velvet Board to have it rolled into a small roll. Now, turn Handle (6) two rounds so that the Roller will deliver a further 2 m.m. of fibres. Draw out the fibres again with Fibre Clip No. 2 and roll them in like manner into a small roll on the Fibre Velvet Board.


Every batch of fibres is then drawn in turn, until all the fibres are clipped up completely. Usually, the fibres are drawn twice for each batch of fibres. When the pointer is set on the Revolvable Dial at degree 16, the Sliding Plate must be lifted up to its original position. Thereafter, when drawing the fibres, the Clip must touch firmly the side close to the edge of the Sliding Plate.

Now the whole cotton tuft is clipped up completely and divided into many batches of fibres. The result is that the lengths of each two neighboring batches are 2 m.m. in difference. Then weigh all batches of fibres separately on a Torsion Balance up to an accuracy of 0.1 milligram to obtain the accurate weight of each batch.

Cotton Sorter—Front View

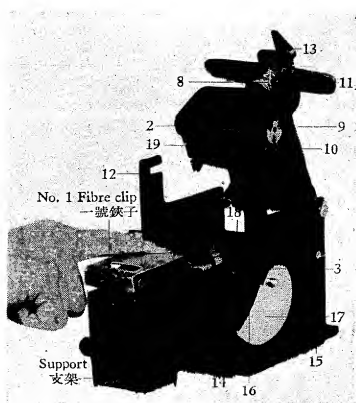Cotton Sorter—Side View

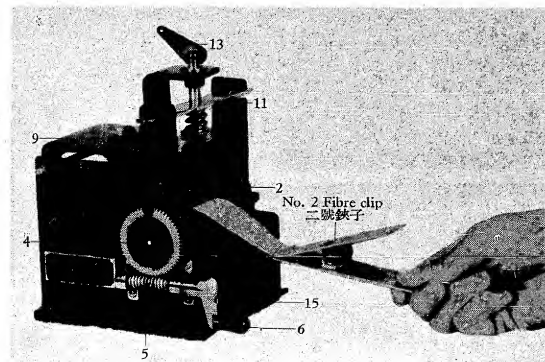Specifications of Drawing Apparatus and Cotton Sorter :—

- | | Drawing Apparatus | Cotton Sorter |
|---------------------------------------------------------------------------------------------------|-------------------|---------------|
| 1. Length: | 155 cm. | 150 cm. |
| Width: | 150 cm. | 120 cm. |
| Height: | 160 cm. | 190 cm. |
| Weight: | 2.3 kg. | 2 kg. |
| 2. Shipping Weight of Drawing Apparatus and Cotton Sorter, including Accessories and Wooden Case: | | 8.3 kg. |
| 3. Hand-operated type. | | |

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMPOT"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.I.).
Cable Address: "INSTRIMP"

TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"

CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

BOOKLET S1006


TEXTILE
TESTING
INSTRUMENTS

STAT

EXPORTED BY:

中國儀器進口公司
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

SHIRLEY ANALYSER

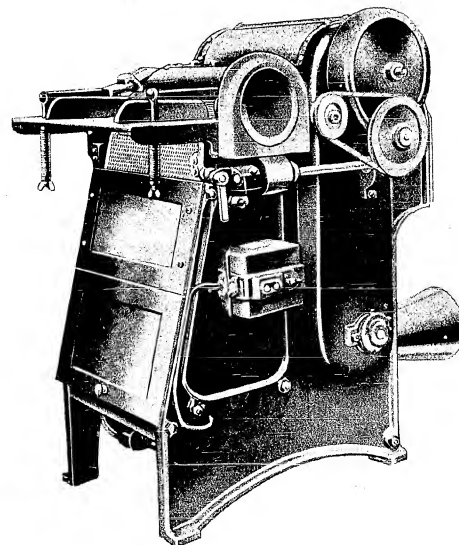Application:-


1. To ascertain percentages of clean cotton and trash of raw cotton or waste.

2. To furnish accurate idea of capability of existing machinery on any particular class or mark of cotton or waste.
3. To determine state of cleanliness of product at any stage in opening, scutching, and carding processes.
4. To ascertain quantity of fibre suitable to be spun in waste discharged by various machines.
5. To determine loss of good fibre in the sequence of opening and scutching.

Overall Dimension

Length 45 inches.
Width 40 inches.
Height 49½ inches.
Weight..... 1433 lbs.

MOTOR DRIVEN LAP METER

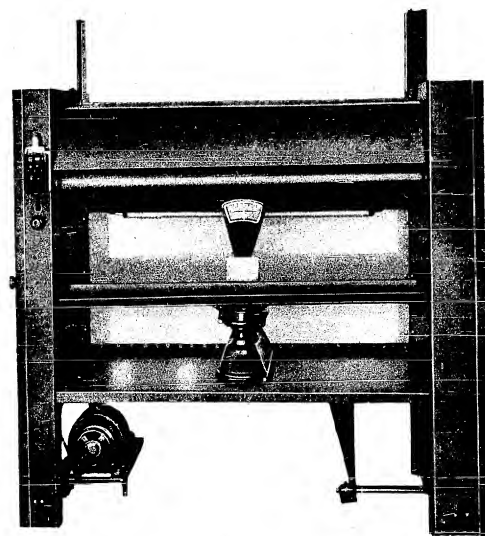Specification:-

1. Testing length - 1 yard.
2. Working width - 40" (45").
3. 1-yard Automatic Stop Motion and Lap Cut-off Device.
4. Attached 1 set 0-26.75 oz. special weighing scale.

5. Wooden roller speed - 2 rpm.
6. Attached one $\frac{1}{4}$ H.P., 380 V, 3-phase motor and one switch.
7. One 220 V fluorescent lamp and one switch.


Direction:

Place lap on iron supporter, start motor and switch on fluorescent lamp, feed lap between wooden roller and $\phi 1\frac{1}{2}$ " roller, and step on the treadle to enable the carrier wheel to engage with the gear on the shaft end of the wooden roller. Having been calendered by the two pairs of rollers, the lap hangs down and its evenness is clearly visible as it passes in front of the light of the fluorescent lamp. As soon as the wooden roller completes one revolution (circumference of which being one yard), the pin of the wooden roller gear pushes away the carrier wheel and so the wooden roller stops. By this time, the lap is cut off by the two pairs of rollers and drops onto the balance dish where it is weighed. Uniformity of weight of every yard is thus ascertained. If yardage is found incorrect, position of $\phi 1\frac{1}{2}$ " roller is to be adjusted.

Overall Dimension

Length	52.4 inches.
Width	26.4 inches.
Height	64 inches.
Weight	1322.5 lbs.

ROVING & SLIVER TESTER (Motor Driven)

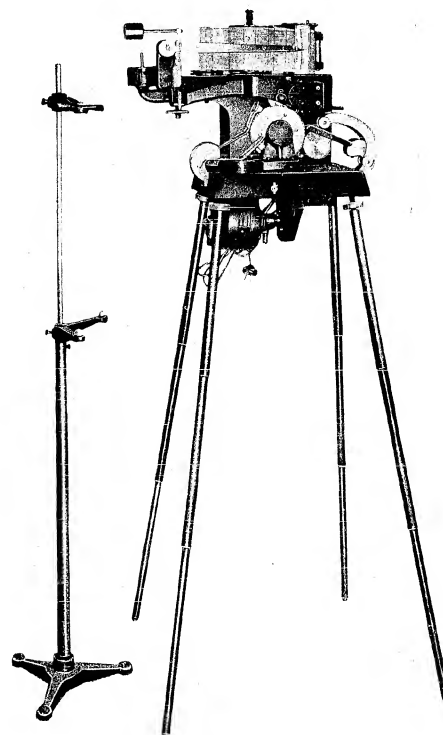
Specification:-

1. Suitable for testing evenness of roving and sliver.
2. Equipped with one 3-phase, 380 or 220 V, 1/8 H.P., 50-cycle, 1420 rpm motor.
3. Attached 6 rolls recording paper and 1 pen.
4. 1 roving stand.
5. Weights: 2½ lbs, for roving, add 1½ lbs. for sliver.
6. Trumpet: For roving 1/16" × 1/2", for sliver 3/32" × 1/2".
7. Forged steel rollers.
8. Space between large and small rollers is detectable under torch light.
9. Space cocupied - 66 × 71 cm.
10. Height - 140 cm.
11. Packing measurement - 64cm × 43cm × 79cm.
12. Gross weight - 86 kgs.
13. Suitable for cotton textile mills.
14. English system.
15. Polished and lacquered in grey.

Overall Dimension

Length	26.4 inches.
Width	28.5 inches.
Height	56 inches.
Weight	189.6 lbs.

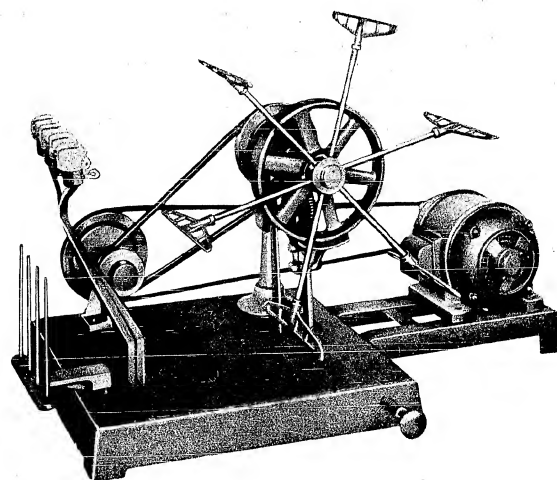MOTOR DRIVEN WRAP REEL**Specification:-**

1. Swift circumference - 1½ yards (or 1 metre for that of metric system).
2. Distance between spindles 2", 4 spindles.
3. Thread guide - ¼" - 1½".
4. Automatic stop motion.

5. Swift speed - 160 rpm.

6. Equipped with one $\frac{1}{4}$ H.P., 220/380 V, 3-phase motor.
1420 rpm.


Direction:

Put bobbins on spindles. Yarn is conducted through thread guide and porcelin cup and is eventually fixed to the hook of one of the brass blades on swift arms. Push knob, traverse rope pulley will connect with swift by means of clutch, and worm on shaft end of swift then drives yardage wheel and swift guide cam. When swift reaches its 80th revolution (100th revolution for the metric system), pin on yardage wheel pushes off the brake hook, enabling clutch to get loosened and the spring to pull the brake lever, and thus the swift is braked. Break yarn at the hook of brass blade marked with colour and remove leas off swift individually.

Overall Dimension

Length	20 inches.
Width	20 inches.
Height	35 inches.
Weight	123.5 lbs.

LEA YARN STRENGTH TESTER (Hydraulic Type)

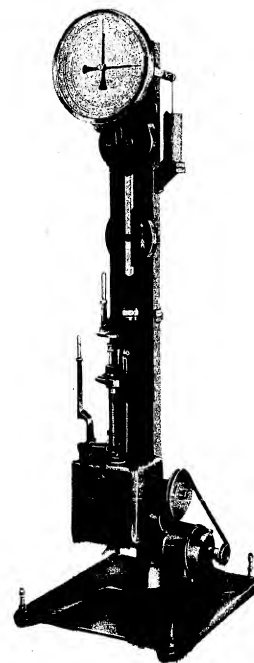Specification:-

1. System-metric or English.
2. Strength range:-
Metric system: Inner graduation - 0 - 100 kgs.
Outer graduation - 0 - 200 kgs.


English system: Inner graduation - 0 - 150 lbs.
Outer graduation - 0 - 300 lbs.

3. Descending speed - 60 cm/min.
4. Swift - 100 cm.
5. Max. descending distance of clamp - 165 cm.
6. Clamp-made of cast iron.
7. Tolerance:-
Inner graduation: 0 - 25 kgs $\pm \frac{1}{4}$ kg, 25 - 65 kgs $\pm \frac{1}{4}$ kg,
65 - 100 kgs $\pm \frac{1}{4}$ kg.
Outer graduation: 0 - 50 kgs $\pm \frac{1}{4}$ kg, 50 - 130 kgs $\pm \frac{1}{4}$ kg,
130 - 200 kgs $\pm \frac{1}{4}$ kg.
8. Levelling position-on sides and rear of U iron.
9. Attached one 3-phase, 380 or 220 V, 50 cycle, 1420 Revolution 0.3 Amp. motor.
10. Attached one V belt.
11. Rotational direction-per direction of arrow on rope pulley.
12. Space occupied - 61cm \times 71cm.
13. Height - 190 cm.
14. Switch-not included.
15. Packing-wooden case, 208cm \times 76cm \times 81cm.
16. Polished and lacquered in grey.
17. Gross weight - 254 kgs.

Overall Dimension

Length 24.4 inches.
Width 28.4 inches.
Height 76 inches.
Weight 560 lbs.

CLOTH STRENGTH TESTER (Hydraulic Type)

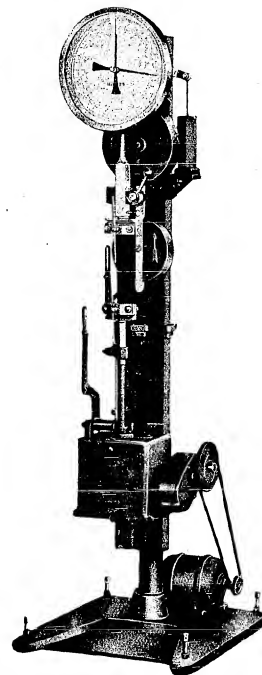Specification:-


1. System-metric or English.
2. Strength range:-
Metric system: Inner graduation - 0 - 100 kgs.
Outer graduation - 0 - 200 kgs.

- English system: Inner graduation - 0 - 150 lbs.
Outer graduation - 0 - 300 lbs.
3. Descending speed - 10 cm/min.
 4. Test cloth width - 5 cm.
 5. Distance between clamps - 20 cm.
 6. Max. descending distance of clamp - 13 cm.
 7. Clamp-made of cast iron.
 8. Tolerance:-
Inner graduation: 0 - 25 kgs $\pm \frac{1}{4}$ kg, 25 - 65 kgs $\pm \frac{1}{4}$ kg,
65 - 100 kgs $\pm \frac{1}{4}$ kg.
Outer graduation: 0 - 50 kgs $\pm \frac{1}{4}$ kg, 50 - 130 kgs $\pm \frac{1}{4}$ kg,
130 - 200 kgs $\pm \frac{1}{4}$ kg.
 9. Levelling position-at sides and rear of U iron.
 10. Attached one 3-phase, 380 or 220 V, 50 cycle, 1420 rpm, 0.3 Amp. motor.
 11. Attached one V belt.
 12. Rotational direction-per direction of arrow on rope pulley.
 13. Space occupied - 61cm \times 71cm.
 14. Height - 190 cm.
 15. Switch-not included.
 16. Packing-wooden case, 208cm \times 76cm \times 81cm.
 17. Polished and lacquered in grey.
 18. Gross weight - 254 kgs.

Overall Dimension

Length 24.4 inches.
Width 28.4 inches.
Height 76 inches.
Weight 560 lbs.

TWIST TESTER

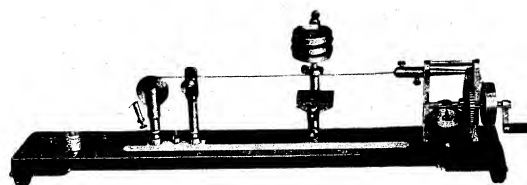
This apparatus is used for determining the twist of cotton, woolen, silk and linen yarns of specific length, fitted with magnifying glass and black and white reflecting plates for showing condition of the yarn while untwisting.

Test length 0-30 cm.
Twist dial 0-500 revolutions (for left and right twists).
Tension weights 5-50 grammes (5 pcs.).
Accuracy 0.01 twist.

Overall Dimension

Length 24.7 inches.
Width 6.5 inches.
Height 7.2 inches.
Weight 14.35 lbs.

COTTON FIBRE TENSILE STRENGTH TESTER

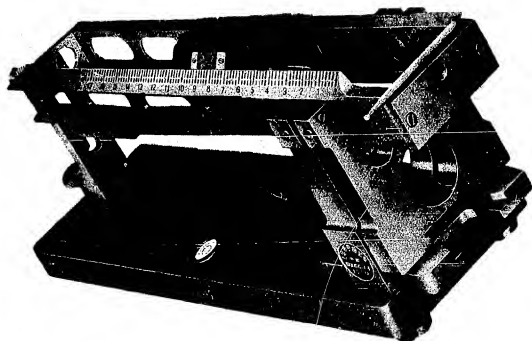This tester is used to determine the tensile strength of a tuft of cotton fibre. Its construction is similar to that of a weighing scale. Being fitted with steel balls, the weight slides easily along the steel rail of the scale lever. When the fibres break, the weight slides to the left and stops on the inclined lever. The point on the scale where the weight stops, indicates the break strength of the fibres.

Capacity	Break strength 10 - 15 lbs.
Testing length	0.464 inch
Accuracy01 lb


Overall dimension:-

Length	15 inches
Width	3.95 inches
Height	6.00 inches
Size of wooden case	L. 19 x W. 7.15 x H. 8.35 inches
Packing measurement	L. 20 x W. 9.2 x H. 10 inches
Net weight	11 lbs
Gross weight	20 lbs

Accessories:-

- 2 - Fibre clamps
- 1 - Sorting device
- 1 - Steel needle comb
- 1 - Socket
- 1 - Knife
- 2 - Black plush boards.

MOTOR DRIVEN TENSION FEELER TYPE TWIST COUNTER

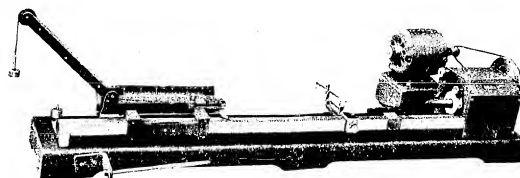This apparatus is adapted to variable testing lengths of 0 - 20 inches. It is provided with a small electric motor and a counter showing from 0 - 9999 turns for left and right. It facilitates the determination of the number of turns of single yarns by opening and twisting again the thread in the opposite direction until the original number of turns is reached as indicated by a small and sensitive tension feeler. Divide the indication of the counter by 2, and the quotient is the number of turns of the tested thread.

Overall dimension (L x W x H) 37" x 14" x 9 1/4"

Net weight 36 1/2 lbs.

Boxed for export 38" x 7 1/2" x 10"

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMPOT"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E.I).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

MICROSCOPES

STAT


CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
SHANGHAI BRANCH
IMPORTERS & EXPORTERS

MODEL "05" BIOLOGICAL MICROSCOPE

A ruggedly constructed instrument designed for general work in the laboratories of educational and research institutes, hospitals, clinics and industries.

This Model "05" Biological Microscope has the following mechanical features:

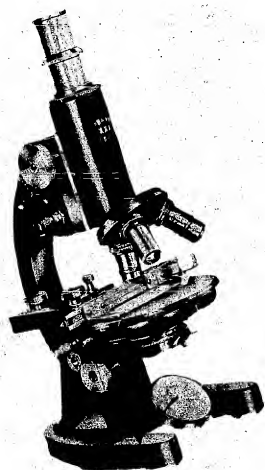Body: Vertical, monocular.
Tube length: 170 mm.

Stage: Round, 115 mm in diameter, with detachable mechanical stage fitted directly into the stage clip holes.

Substage: Rack and pinion with Abbe Condenser Iris diaphragm.

Focusing Adjustment: Standard rack and pinion coarse adjustment. Standard side lever type fine adjustment, 0.01 mm.

Nose Piece: Revolving, dust-proof, centered and parfocalized.

Finish: In rich, lustrous black. Exposed metal parts are chromium plated.

Optical Equipment:

Huygenian Eyepieces:	Achromatic Objectives:	Condenser:
5x 10x 15x	8x 0.25 N.A. 45x 0.63 N.A. 100x 1.32 N.A. (Oil immersion)	Abbe Condenser: 1.30 N.A.

Accessories:

Spanner	1 piece	Cedar Wood Oil	1 bottle
Stage Spring Clips	1 pair	Filter Glass	1 piece


Packed in a handsome hardwood case of 200 x 240 x 360 mm.

Gross weight: 9.20 kgs. 上海儀器廠製造

Model "08" RESEARCH MICROSCOPE.

Best suited for research institutes, clinics, hospitals and universities.

This Model "08" Research Microscope has the following mechanical features:

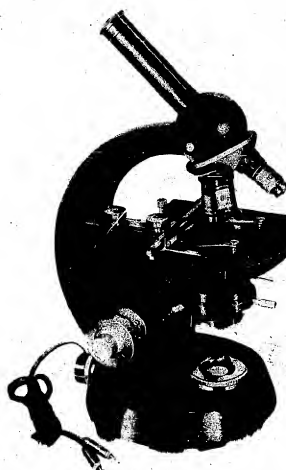Body: Inclined Monocular, inclination 45°, tube length 160 mm.

Stage: Fixed round mechanical stage.

Substage: Rack and pinion with Condenser (N.A. 1.30). Iris diaphragm.

Focusing Adjustment: Coarse adjustment by means of rack and pinion. Fine adjustment by means of ball-bearing device, 0.002 mm graduation.

Nose-piece: Revolving, quadruple.

Objectives: 8x, 45x and 100x (oil immersion), equipped with built-in spring device to avoid damage that may otherwise occur to objectives and microscopic slides or specimens when contacting inadvertently.


Illumination: Built-in Electric illuminating device.

Optical Equipment:

Achromatic Objectives:		Huygens Eye-Pieces:	Total Magnification.
Dry:	Oil Immersion:		
8x 45x	100x (N.A. 1.32)	5x, 10x, 15x.	40x — 1500x.

Packed in handsome varnished cabinet.

Gross Weight: 13.5 kgs.

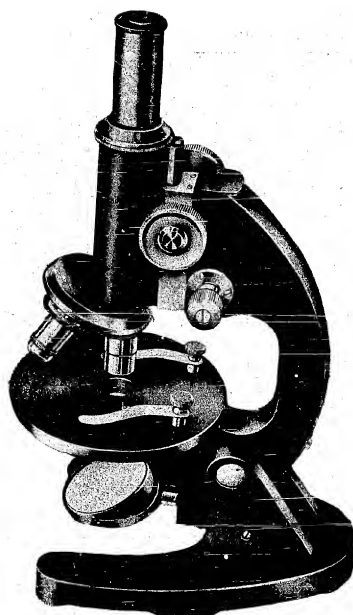Model "01"
BIOLOGICAL MICROSCOPE

Best suited for school experiment and laboratory purpose.

Model "01" Biological Microscope has the following mechanical features:

Body: Vertical, monocular. Tube length 170 mm.

Stage: Round 120 in diameter, with revolvable quintuple diaphragm.

Focusing Adjustment: Coarse and fine adjustments.

Nose-piece: Revolving, triple.

Optical Equipment:

Achromatic Objectives	Huygens Eye-pieces	Total Magnification
8x, 45x	5x, 10x, 15x	40x - 675x

Packed in handsome varnished cabinet.

Gross weight: 5 Kgs.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:
Rr-Li-Kou, 11st China, Peking, China.
Cable Address: "INSTRIMP"
AGENTS:
HONGKONG: China Resources Co.,
Bank of China building 11/12 floor, Hongkong.
Cable Address: "CHIRCO HONGKONG"

BRANCH OFFICES:
SHANGHAI: 27 Chongshin Road (F. 1).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

TORSION BALANCE

STAT

中國儀器進口公司
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(INCORPORATED IN HONGKONG)

TORSION BALANCE


Model TW-TB

SPECIAL FEATURES:

- Highest accuracy, sensitivity and reliability.
- No bearing friction.
- With clear and readable scale.
- Anti-parallax reading.
- Rapid weighing.
- Convenient adjustment.
- Insensitive to external influences of magnetic forces.

The Torsion Balance is supplied in the following capacities:

Measuring Range in mgm.	Each division in mgm.	Reading in mgm.	Cat. No.
0-5	0.01	0.005	STW 1001
0-10	0.025	0.01	STW 1002
0-25	0.05	0.025	STW 1003
0-50	0.1	0.05	STW 1004
0-100	0.25	0.1	STW 1005
0-250	0.5	0.25	STW 1006
0-500	1.0	0.5	STW 1007
0-1000	2.0	1.0	STW 1008

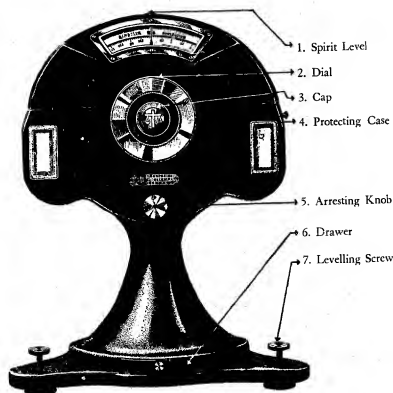Torsion Balance Model TW-TB

METHOD OF OPERATION

INSTRUCTION FOR SETTING UP:

1. Put the balance in a definite firm place where the weighing is done.
2. Level the balance by adjusting the leveling screw (7) until the bubble in the spirit level (1) is centered.
3. Open the protecting case (4) and cautiously hook up the weighing pan. Then close the case.
4. Set the zero point of the dial on the hair line by turning the dial (2). Turn the arresting knob (5) to the left for about 90° to release the balance and check if the hair line, the pointer and the reflection of the pointer in the mirror all fall in a line. When this is so, arrest the balance by means of turning the arresting knob to the right for about 90°.

Should the hair line, the pointer and the reflection of the pointer in the mirror do not fall in one line at this stage, open the cap (3) and, by using the supplied screw driver, the adjustment is made by turning the screw, which is located under the cap (3), to the right if the pointer falls off to the left of the hair line, or vice versa, until the pointer is in the proper position.

5. Open the protecting case and place the check weight on the weighing pan by means of the supplied tweezers, then close the protecting case again. Set the full mark of the scale on the hair line, then release the balance by turning the arrester to the left to check if the hair line, the pointer and the reflection of the pointer in the mirror all fall again in one line.
6. Arrest the balance again by turning the arresting knob to the right. Open the protecting case and take away the check weight from the weighing pan. Finally, set the zero point of the dial on the hair line.

NOW, THE BALANCE IS READY FOR USE.

WHEN IN USE:

1. Cautiously place the weighing material on the weighing pan and close the protecting case to shield the hook and the sample from draught. First turn the arresting knob to the left to release the balance, then turn the dial until the pointer falls in line with the hair line and the reflection of the pointer in the mirror, whereby the accurate weight is directly indicated in the scale by the pointer, and the weighing is completed.

Standard Equipment: (kept in drawer)

- 1 screw driver.
- 1 pair of tweezers.
- 1 checking weight for the capacity of the balance.

FIELDS OF USES FOR THE TORSION BALANCE:

The Torsion Balance is widely used in Research Laboratories, Educational Institutes and different lines of industries. It serves:

1. To check the weight-alterations due to evaporation and incineration processes.
2. To examine seeds of all kinds.
3. To investigate raw fibres, yarn and fabric samples, raw materials and finished products.
4. To determine the hygroscopy and moisture contents.
5. To determine the weight of:
 - paper samples, medicine powder, tablets, pills, incandecent filaments, etc.

Dimensions of carrying case 32.5 x 24 x 40.5 cm.
 Gross Weight (including carrying case) 8 kg.
 Net Weight 5 kg.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION


HEAD OFFICE:
 Er-Li-Kou, Hsi Chiao, Peking, China.
 Cable Address: "INSTRIMPOT"

BRANCH OFFICES:
 SHANGHAI: 27 Chungshan Road (E.I.).
 Cable Address: "INSTRIMP"
 TIENTSIN: 171 Kien Shieh Road.
 Cable Address: "INSTRIMP"
 CANTON: 25 Tai Ping Road, S.
 Cable Address: "INSTRIMP"

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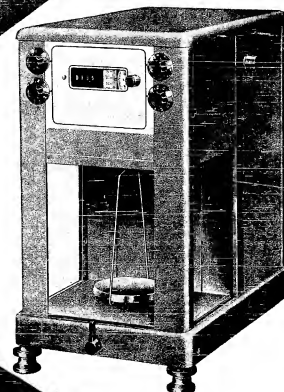Analytical Balances

Representative products of our modern industries, enjoying wide admiration for their high sensitivity & durability

Exported by:

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

SINGLE-PAN PROJECTION READING ANALYTICAL BALANCE

Latest scientific researches have enabled us to offer to the public an entirely new model of analytical balance, the construction of which is made possible by the employment of a principle vastly different from that hitherto used. Instead of the conventional 3 knife edges and 2 pans, the new balance has only 2 fine agate knife edges and 1 pan. Its weights are loaded on a support suspended to one end of the beam with a counter weight fixed to the other end. These weights can be readily detached from or reloaded onto the support by means of remote control knobs constructed outside the case enclosing the balance. Weighing is done by simple manipulation of the knobs by detaching the equivalent weights from the support and with the help of the electric projected vernier, the accuracy of reading therefrom can reach 0.05 mg. The result of the weighing can be easily readable from the optical scale specially provided for this purpose. The capacity of this balance is 100 grams.

The balance is equipped with a built-in transformer adaptable to 110v, 130v, 160v, 220v A.C. and is enclosed in a sturdy aluminium alloy case protected with anti-oxidizing laquer, with glass sliding doors and an overall size of 440×410×270 mm, weighing about 17.5 kilograms.

ADVANTAGES:

1. Accuracy:

- a. All weights are fixed inside the balance and handled with ease by control knobs outside the case. They are chrome plated to resist corrosion and thus unalterable.
- b. Constant loading on the beam eliminates bending errors and ensures uniform sensitivity up to the maximum loading.
- c. The optical reading range is 115 milligrams. All weighing under 100 milligrams is absolutely accurate.
- d. The weighing indicator scale is fixed on the front control panel.


(figure 1) and it is divided into three parts:

From 0.1 gram to 99.9 grams read directly from indicator (1)
 From 1 milligram to 100 milligrams read directly from projection scale (2)
 From 0.05 to 0.1 milligram read directly from vernier (3)

2. Rapid Weighing:

The time required for weighing with this balance is less than 1/3rd of that of an ordinary balance. After the zero point is adjusted, weighing is done simply by turning the control knobs. Every weighing can be done in less than 30 seconds. Due to the ingeniously designed air-damping device, the balance is brought to rest quickly after only a few oscillations.

Model NUCSP

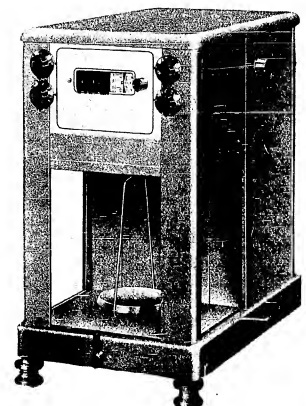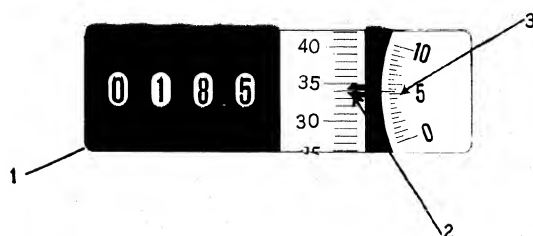
Figure 1



3. Simplicity:

- a. No arithmetic calculation is needed since the total weight is read directly off the indicator scale.
- b. When weighing, simply detach weights corresponding to the article weighed by turning the control knobs, and the total weight is read right off the projection scale.

DIAGRAM:

1. Adjustable plates for the elimination of errors due to atmospheric changes
2. Air damping device
3. Arrestment lever
4. Beam arrestment
5. Center knife edge and plane of agate
6. Compensating weight for pan and weights
7. Compensating weight for adjustment of the zero point
8. Fine agate knife edge and plane
9. Optical scale plate
10. Pan brake lever
11. Pan, circumference 100 mm.
12. Sensitivity compensation adjustment weight
13. Set of weights; two 25-g; two 10-g; five 5-g; one 2-g; two 1-g; one 500-mg; one 200-mg; and two 100-mg.
14. Stirrup from which the weights and pan are suspended

Diagram:

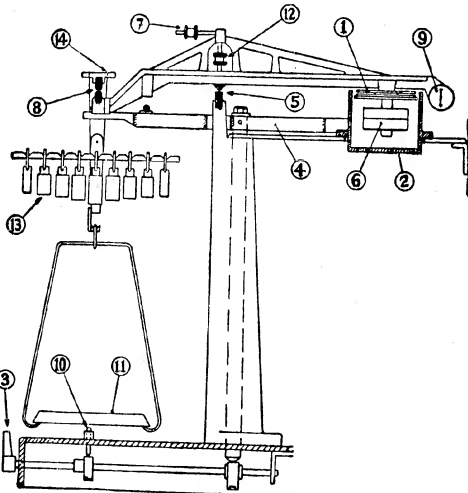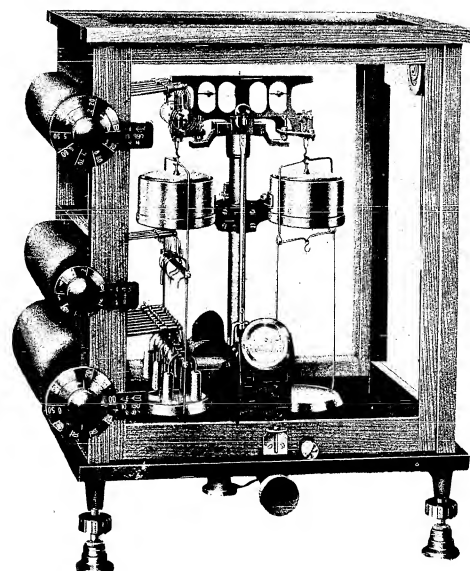**CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION**

Model NCAPD

ANALYTICAL BALANCE
AUTOMATIC PROJECTION READING FOR RAPID AND RELIABLE WEIGHING

Capacity: 200 grams

Sensitivity: 0.1 milligram

SPECIAL FEATURES:

- Air-damping
- Projection Reading
- External Application of All Weights
- Zero Point Adjustment from Outside of Case

CONSTRUCTIONS:

Beam:-

Brass alloy in flat form, 140 mm long, chrome plated

Knife edges and plane:-

Finest agate

Weights:-

Brass, chrome plated

Pans:-

Slightly concave, chrome plated

Compensating suspension:-

Movable in all directions, pendulum motions of the pans without any effect upon the blocks

Air-damping:-

Without diminution of sensitivity, brings the balance to rest in 3 to 4 seconds

Projection reading with automatic light switch and adjustable device:-

Range 0.1 mg. to 10 mg, direct reading to the exactness of 0.1 mg. A.C. transformer included, applicable to house current

Mechanical manipulation of weights from outside the case:-

From 10 to 990 mg. by means of the upper dial &/or turning knob

From 1—9 g. by means of the middle dial

From 10 to 190 g. by means of the lowest dial &/or turning knob

Arrestment:

By means of central knob which also operates as switch for projection light


Zero point adjustment:

From outside the case by means of a small knob beside the arrestment knob

Base plate:

Mirror glass with levelling screw and shock absorbers

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

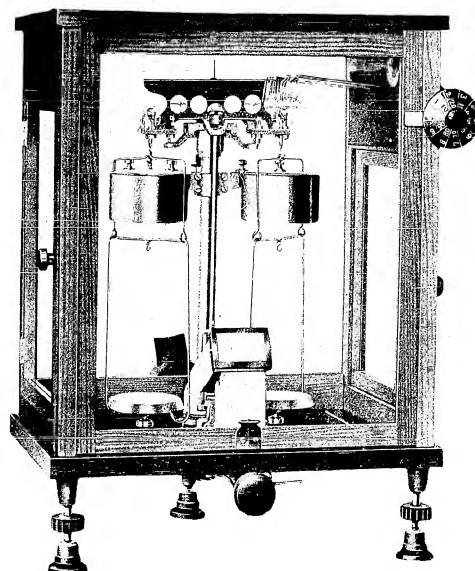Model SA00G


ANALYTICAL BALANCE **SEMI-AUTOMATIC PROJECTION READING FOR RAPID & RELIABLE WEIGHING**

Capacity: 200 grams

Sensitivity: 0.1 milligram

With all special features and similar constructions as NCAPD except that the mechanical manipulation of weights from outside the case is only for range of 10 mg. to 990 mg.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

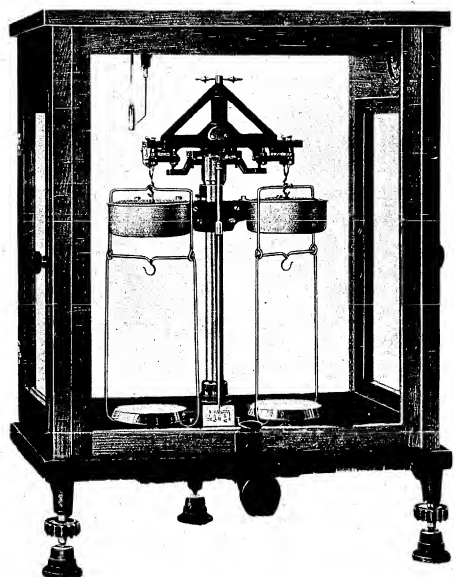Model FANOF

ANALYTICAL BALANCE

Capacity: 100/200 grams

Sensitivity: 0.4 milligram

Beam:

140 mm long, aluminium alloy with colorless anti-oxidizing lacquer

Graduation:

To ensure accurate placement of rider, the indicating scale attaching to the beam is serrated and notched into 50 divisions from Central "0" to each side, each division representing 0.2 mg when placing a 10 mg rider on the scale

Knife edge & plane:

Agate

Rider slide:

Vertical type, spring loaded, traversing the whole length of the beam

Releasing Mechanism:

Automatically acting on hangers, beam and pans, operated from outside the case

Pans:

Chrome plated, slightly concave, diameter 73 mm

Pan supporters:

Adjustable with special spring plunger so that oscillation does not occur when setting the balance into operation

Reading Index:

20 divisions from "0" on the right to "20" on the left

Air-damping:

The swing of beam is damped by detachable air damping device


Finish:

All metal parts are protected with either chrome plating or anti-corrosion lacquer

Case:

Polished teak, with counterpoised iron slide and hinged side doors. Black plate-glass or marble base. Levelling screw and bubble levelling system

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

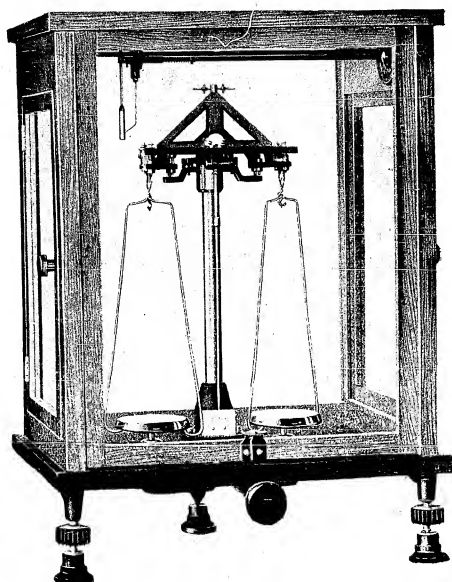
Model NODOF

ANALYTICAL BALANCE

Capacity: 100/200 grams Sensitivity: 0.4 milligram

This balance is similarly constructed as FANOF with the exception that no Air-damping device being provided.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

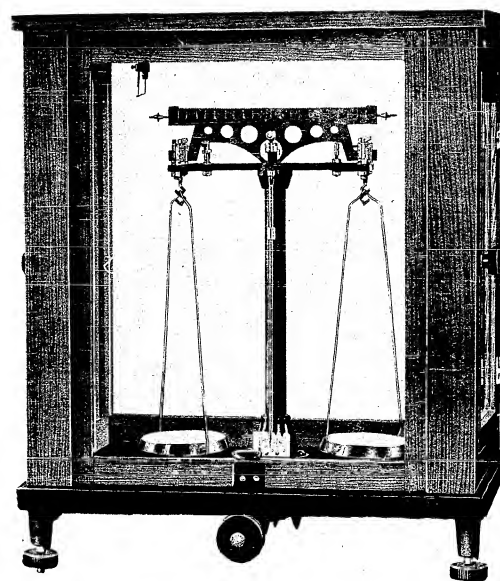
Model SNBOM

PRECISION BALANCE

Capacity: 100/200 grams Sensitivity: 1 milligram

A precision balance for schools, laboratories and/or chemical works. With rider slide. Agate knife edges and planes. Arrestment controlled by central knob.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

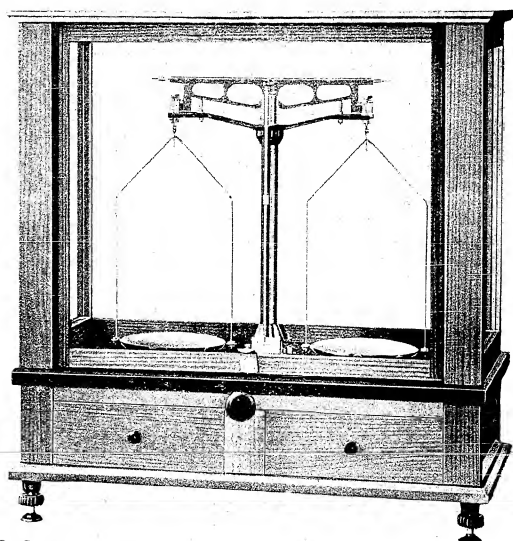
Model CTNAO Sensitivity: 4 mg.
Model CTFMG Sensitivity: 20 mg.

PRECISION BALANCE

Capacity: 100 grams. Sensitivity: 4 mg. or 20 mg.

This precision balance is specially designed for textile industries, with brass alloy beam, riderless, 2 concave bakelite round pans. These balances are contained in polished teak cases with counterpoised front slide doors, black lacquered wooden base, levelling screw and plummet.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

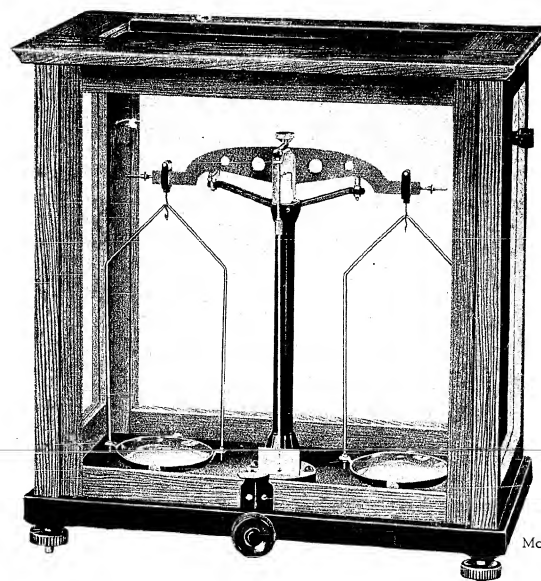
Model MAWTI

PRECISION BALANCE

Capacity: 100 grams Sensitivity: 10 milligrams

This precision balance is specially designed for mineral industries, schools and laboratories, riderless, concave bakelite pans diameter about 75 mm, aluminium alloy beam with colorless anti-oxidizing lacquer, agate knife edges and planes. The case is of polished hardwood, counterpoised front slide door, black lacquered wooden base. Levelling screw and plummet.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

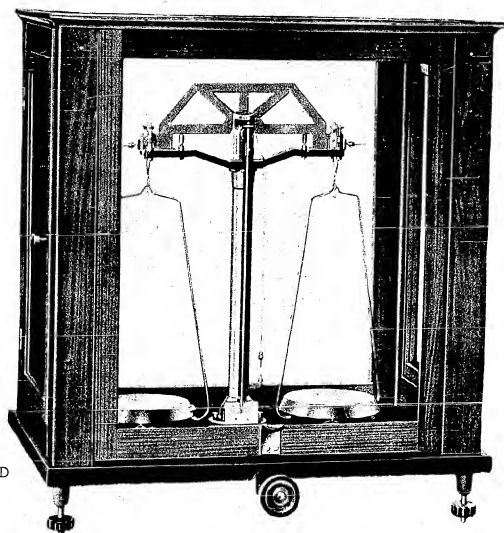Model IABHD

INDUSTRIAL ANALYTICAL BALANCE-HEAVY DUTY

Capacity 1000 grams


Sensitivity 2 milligrams

Diameter of Pans 107 mm

Height of Pan Holders 255 mm

Length of Beam 270 mm

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

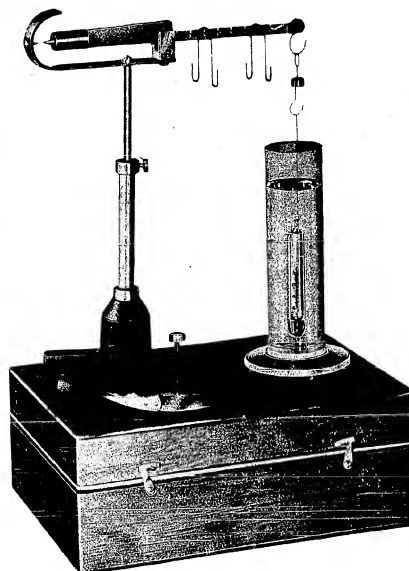Table for DENSITY OF DISTILLED WATER
at different temperature

Temper- ature	Density	Temper- ature	Density
0	0.999,867	21	0.998,017
1	0.999,927	22	0.997,795
2	0.999,968	23	0.997,563
3	0.999,992	24	0.997,321
4	1.000,000	25	0.997,069
5	0.999,992	26	0.996,808
6	0.999,968	27	0.996,538
7	0.999,929	28	0.996,258
8	0.999,876	29	0.995,969
9	0.999,808	30	0.995,672
10	0.999,727	31	0.995,366
11	0.999,632	32	0.995,052
12	0.999,524	33	0.994,728
13	0.999,404	34	0.994,397
14	0.999,271	35	0.994,058
15	0.999,126	36	0.993,711
16	0.998,969	37	0.993,356
17	0.998,801	38	0.992,993
18	0.998,621	39	0.992,622
19	0.998,430	40	0.992,244
20	0.998,229		

Model SGBTF

SPECIFIC GRAVITY BALANCE-WESTPHAL PRINCIPLE

Complete with rider weights, plummet (displacement 5 cc) with platinum wire and stainless steel hook, the thermometer connected therewith is graduated from 0 to 35 degrees Centigrade and it is for determination of Specific Gravity of all liquids up to 2 S.G. The beam of the balance is fitted with steel knife edge supported in agate channel and is precisely serrated in 10 divisions. Before leaving the factory, the balance is subjected to rigid test to conform strictly to the table herewith.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

ANALYTICAL WEIGHTS

Calibrated weights of high accuracy are supplied together with our balances without extra charges. They are placed in velvet lined polished wooden box in the sequences:

100g	1	500 mg	1
50g	1	200 mg	2
20g	2	100 mg	1
10g	1	50 mg	1
5g	1	20 mg	2
2g	2	10 mg	1
1g	1	5 mg	1
		2 mg	2
		1 mg	1

The weights of 50 mg and under are made of aluminium.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

Er-Li-Kou, Hsi China, Peking, China.
Cable Address: "INSTRIMP"

BRANCH OFFICES:

SHANGHAI: 27 Chungshan Road (E. 1).
Cable Address: "INSTRIMP"
TIENTSIN: 171 Kien Shieh Road.
Cable Address: "INSTRIMP"
CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

20

STAT

INSTRUCTIONS FOR WORKING AND ASSEMBLY OF OUR ANALYTICAL BALANCES

中國儀器進口公司
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
SHANGHAI, CHINA IMPORTERS & EXPORTERS

AUTOMATIC PROJECTION READING ANALYTICAL BALANCE (SINGLE PAN TYPE)

DIRECTIONS FOR USE:

"Single Pan" Automatic Projection Reading Analytical Balances are instruments of highest precision. Although they are simple in operation, it is essential they are handled with utmost care. The following rules should therefore be observed:

1. Do not load or unload the pan until the balance is arrested.
2. Move the release slowly until the light goes on.
3. Handle the rider bar gently.
4. When the balance is arrested, rotate the left knob clockwise from 10 to 90 g., the two right knobs counter clockwise from 1 to 9 g. and 0.1 to 0.9 g. respectively.

WEIGHING:

Before each weighing, level the balance by means of the spirit level (28) provided within the case and adjust the zero position as follows:

With no load on the pan, the knobs are set at 000.0 and the micrometer is adjusted to zero. Then move the arrest lever slowly to the left. Let the beam stabilize. By turning knob (24) on the right side of the balance, the zero on the projected scale should be set exactly between the fixed double-line (25) on the screen. If the above adjustment is not sufficient, remove the cover, check the beam (15) and suspension gear (5) to locate the cause of friction, and then move the horizontal adjustment weight (7) on the beam accordingly. Move the weight (7) forward if the projected scale is too high, and backward if it is too low.

When the above adjustment has been completed, a trial weighing should be made as follows:

Place a known weight, say 0.1 g. on pan. Partially release the pan by turning the arrest lever to the right. If the weight shown on the projected scale is lighter than 100 mg., projected scale and pan should be slightly lifted; if heavier, slightly lowered. When the scale reads exactly 0.1 g., move arrest lever slowly to the left to give full release of the beam. A correct result will thus be obtained. Example of weighing:

Suppose the weight of the sample is 18.53452 g. Set knob A at 10 g. Turn arrest lever to the right, the projected scale will move downwards, i.e. its zero remains about 5 divisions below the double-line on the screen. This shows the sample on pan weighs more than 10 g. Arrest the balance, set knob A at 20 g., and again turn arrest lever to the right. The scale moves upwards. The weight of the sample must therefore be between 10 and 20 g. Arrest balance, switch knob

A back to 10 g. Move arrest lever to the right and turn knob B until the scale moves upwards. This will happen at 9 g. At this moment, the counter indicates 18.5 g. The remaining 4 decimals are determined optically by turning the arrest lever fully to the left. After a few small oscillations, the optical scale comes to rest between the 34th and the 35th divisions. The small knob (27) is now turned until the line of the 34th division of the projected scale is set exactly between the double-line on the screen, as in fig. II. The final result will read 18.53452 (see figure III).

Return knob B to 8 g. and proceed likewise with knob C. Result 0.5g.

SENSITIVITY:

To check the sensitivity of the projected scale, place a 100 mg. weight on the pan, turn knob C to 0.1 g., release balance, and adjust zero. (Cover and windows must be closed). Before the projected scale correctly shows 100 mg., do not arrest and do not turn knob from 0.1 g. back to 0. If the scale reads 99 for instance, the sensitivity must be increased by raising the vertical adjustment weight (6). If it reads more than 100, the weight should be lowered accordingly. The balance must have reached the room temperature before the sensitivity can be adjusted.

IMAGE:

If the image of the optical scale is not clear enough, remove the cover, loosen screw (22a) and move objective (22). Tighten screw (22a) again.

If the line of the projected scale do not coincide with the double-line of the screen, rotate prism turntable (23).

If the projected image is coloured or too dark, move bulb (19) upwards.

PAN:

When the balance is arrested, the pan should just touch the brake (1) below it. If the braking is too great or too little, raise or lower the pan by bending or unbending suspender (15).

In case the pan requires to be replaced by one of different weight, screw (25) may be added onto pan-bridges to restore balance.

FOUL CONTACT:

When the balance is released, the projected scale should swing up and down 2 to 3 times before it comes to rest. Coming to rest abruptly indicates that the beam or the suspension is fouling somewhere. Remove the cover to locate cause of friction.

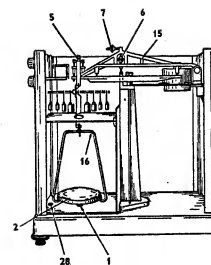
FIGURE I.

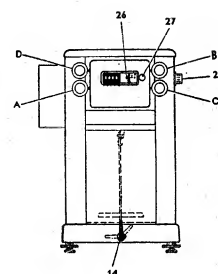
FIGURE II.

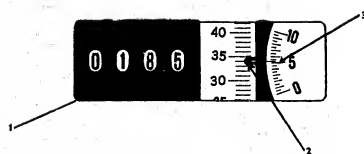
FIGURE III.

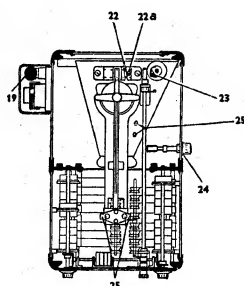
FIGURE IV.

AUTOMATIC ANALYTICAL BALANCE

for

Rapid and Reliable Weighing

Capacity: 200 gm. Sensitivity: 0.1 mg.

INSTRUCTIONS FOR ASSEMBLY:

1. Level the balance with the help of the spirit level by turning the levelling screws.
2. Ascertain the type of current and voltage to be adopted before connecting the transformer.
3. Fix the lamp housing on the condenser tube in the rear, insert one end plug into the transformer and another into the socket underneath the balance base plate.
4. Insert arrest knob into shaft. By turning the knob to right or left, the beam support comes downwards. Grasp pointer of beam at its middle, clean knife edges gently with fine soft cloth or handkerchief, then mount beam on support. Turn the knob to right or left to raise the support.
5. Suspend damping cylinders from the lower hooks of suspension pieces.
6. Suspend pans by inserting into upper hooks of stirrups, after the arrest device has been fitted on. If the arrestment is too great or too little, lower or raise the device by turning its cap.
7. If the image is not bright enough, loosen the knurled-screw which holds the lamp socket in place and shift the socket until the scale image becomes brightly and uniformly illuminated. Then fasten the knurled-screw tightly. If the scale image is not perfectly focussed, loosen set screw of objective and shift same along its two guiding rods until the scale is sufficiently clear.
8. Release the balance, and check whether the pointer coincides with the graticule. If any adjustment is required, it can be effected by means of the screws on the beam.

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INTRODUCTION

1. The material used is an alloy of nickel and bronze guaranteed, to last a lifetime and keep the original good look permanently.
2. All single parts are precisionally made, all fit perfectly.
3. All pen parts are made of high grade steel with scientific heat treat process, assuring to draw lines within 1200 meters without resharpening.
4. All individual instruments have been tested before leaving the works to assure high accuracy.
5. The parts of the screw threads and accessories are made of stainless steel so that they never oxidize and give long service. Others are heavily chromed.
6. The cases from 2001 to 114 are flat and smooth with rounded corners and edges, covered with black fine oil cloth and lined with good quality velvet.
7. The cases of 115 and 121 are made of fine solid wood and lined with good quality velvet, keeping the appearance handsome and durable.


2001

1 : 2

6" compass with pen, pencil, divider and extension bar, all interchangeable.

Precisionally used for drawing circles about 4" in diameter up to 12".


Both legs of the compass are supplied with joint or hinge, so that the points can be perpendicular to the paper.

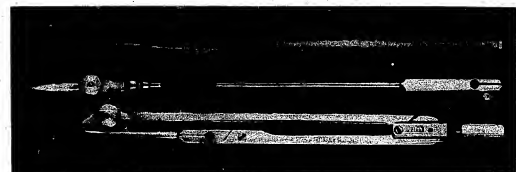100

1 : 2


6" compass with pen, pencil, divider, extension bar and a handle for ruling pen attachment, thus forming a drawing pen. Suitable for beginners and users during travelling.

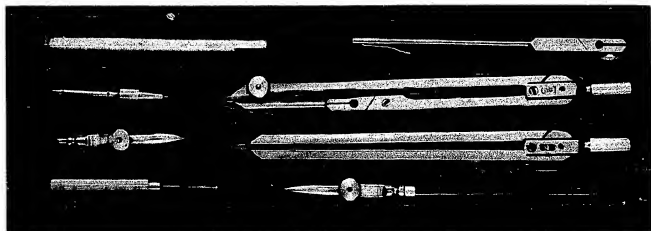103

1:2


School drawing set, consisting of large compass with attachments, divider, ruling pen, handle for pen attachment and screwdriver.

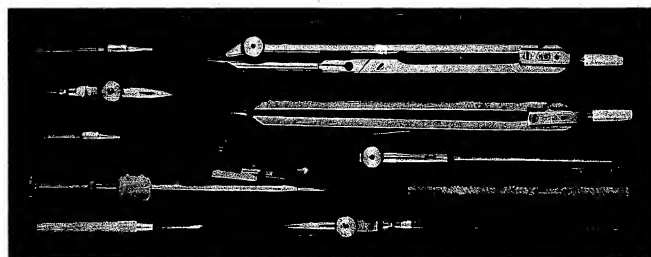104

1:2


Including all the individual instruments in 103, with the addition of a drop compass which is used only for very small circles, having the merit of time-saving and accuracy.

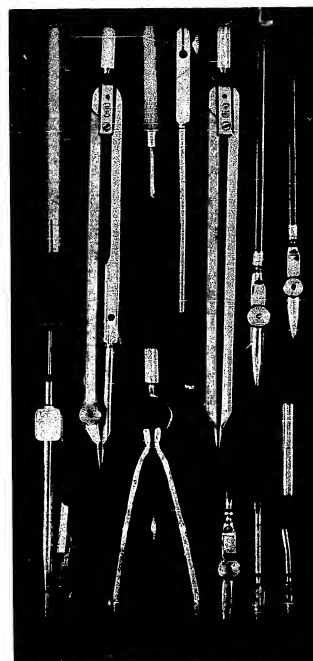106

1:2


Very attractive students' set, consisting of 6" large compass with attachments, divider, spring bow divider, drop compass, small and medium size ruling pens, handle for pen attachment, screwdriver and lead box.

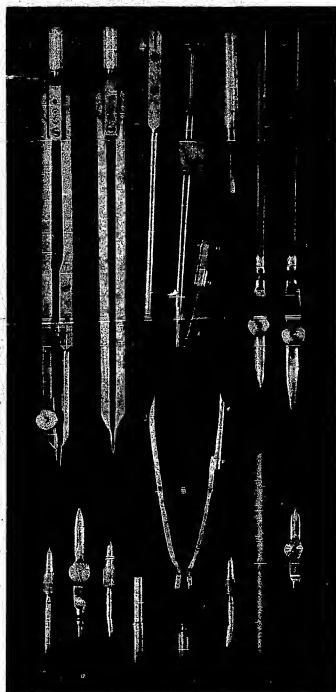107

1:2


Another attractive students' set, consisting of some instruments as 106, but with spring bow attachments instead of spring bow divider and large size ruling pen instead of small size.

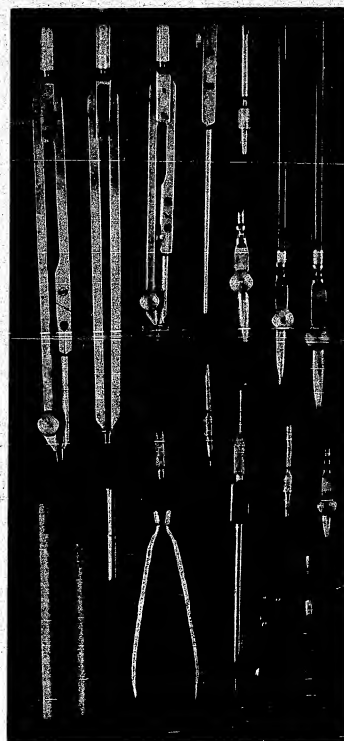109

1:2


Suitable drawing set for engineers and technicians, consisting of two compasses with attachments, two dividers, three ruling pens, one double-spring drop compass, two handles for pen attachments.

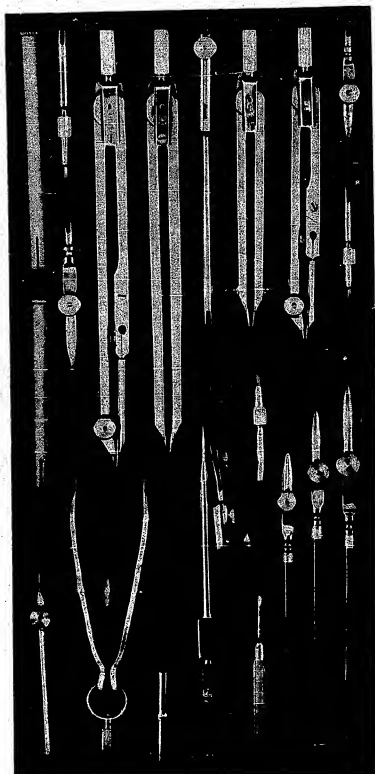111

1:2


A special combination set for engineers and architects, consisting of two compasses with attachments, three ruling pens, one drop compass and accessories. The concentric piece is provided to avoid leaving a large visible hole when drawing a series of concentric circles.

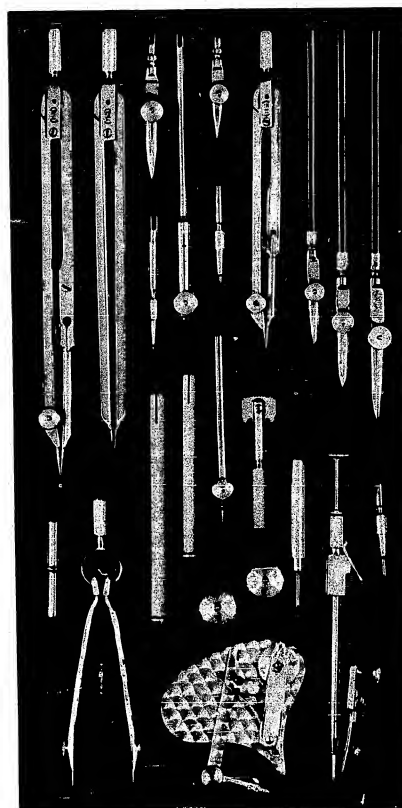114

1:2


A comprehensive precision drawing set for draftsmen, architects, surveyors and engineers, consisting of two compasses, two dividers, three ruling pens, drop compass, dotting instrument and other accessories.

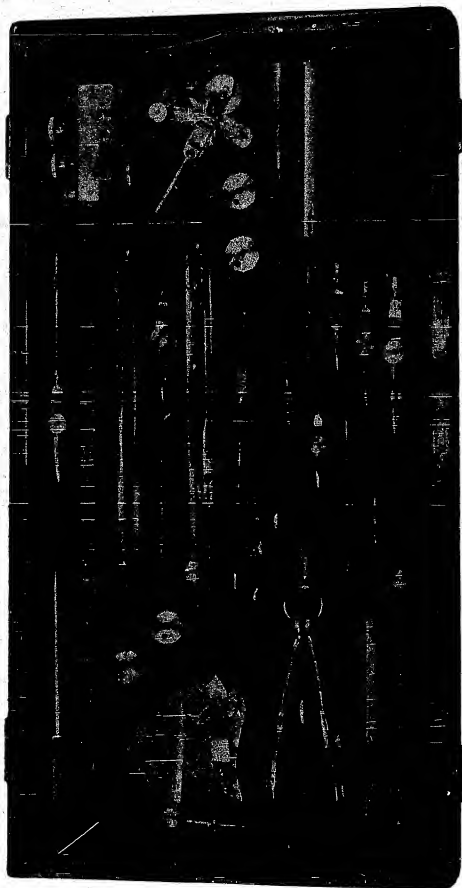115

1:2

Special set for engineers, architects, technicians and draftsmen, consisting of compasses, divider, spring bow compass, drop compass, beam compass, two dotting instruments, three ruling pens and other attachments and accessories.

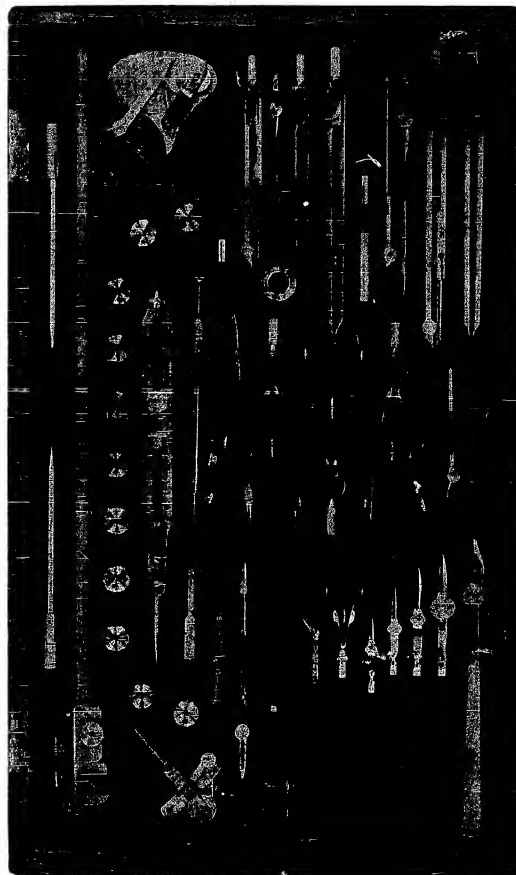
121


1:25

The complete precision drawing set assorted with so many special instruments such as beam compass, dotting instruments, curve pens, border pens, proportional divider and hairspring divider that it is greatly appreciated by the modern engineers, architects, technicians and draftsmen.


SPRING BOW COMPASS

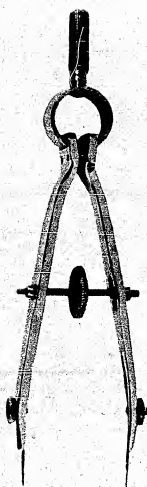
2007

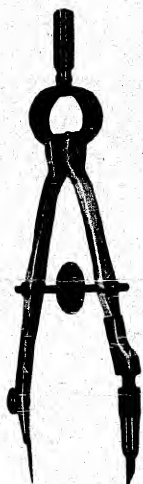2008

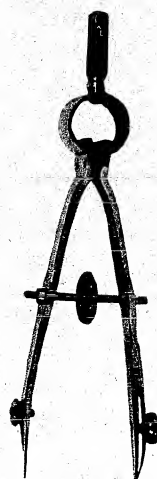
2009

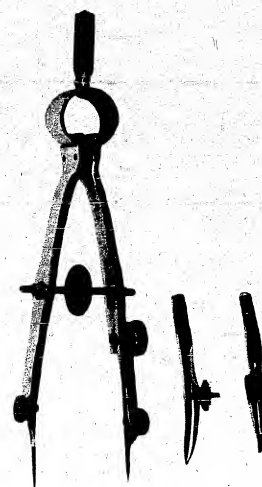2010

The legs of the spring bow compass are fixed to a spring head and regulated by means of a nut between the legs with adjusting screws at both sides. The merit is that the distance can be easily adjusted as required and once adjusted, the legs will never move no matter how many lines are divided or arcs drawn.


SPRING BOW COMPASS

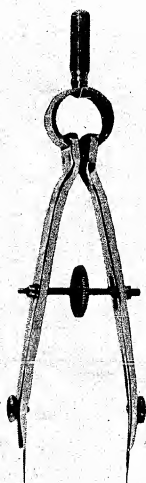
2007

2008

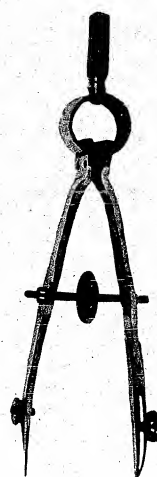
2009

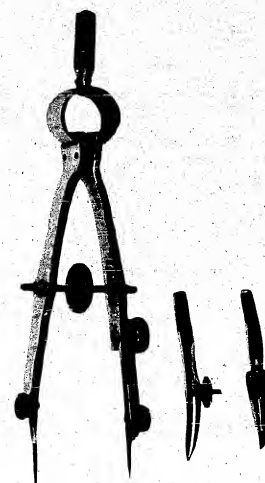2010

The legs of the spring bow compass are fixed to a spring head and regulated by means of a nut between the legs with adjusting screws at both sides. The merit is that the distance can be easily adjusted as required and once adjusted, the legs will never move no matter how many lines are divided or arcs drawn.

COMPASS

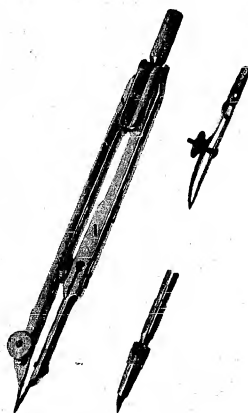
2004

4" compass with pen, pencil and divider attachments. Same construction as 2001. Suitable for circles below 4" in diameter.

DIVIDER

Similar to the compass, with points on both legs, but without joints.

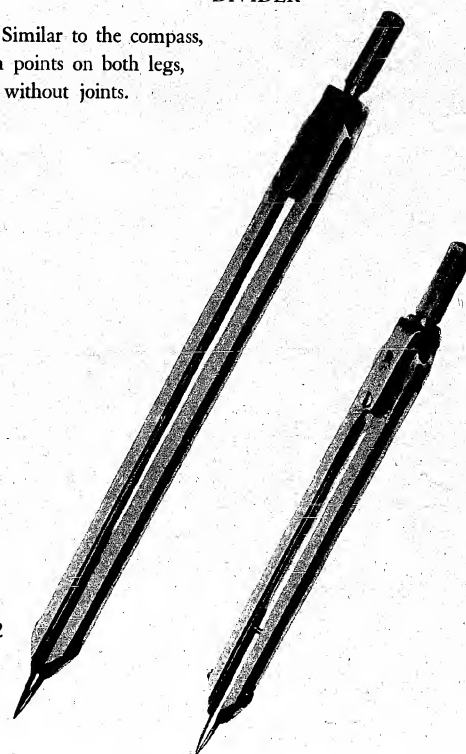
2002

2005

CROSS JOINT RULING PEN

The special advantages of these pens: convenience for cleaning; keeping the lines in same thickness when returned to the original position.

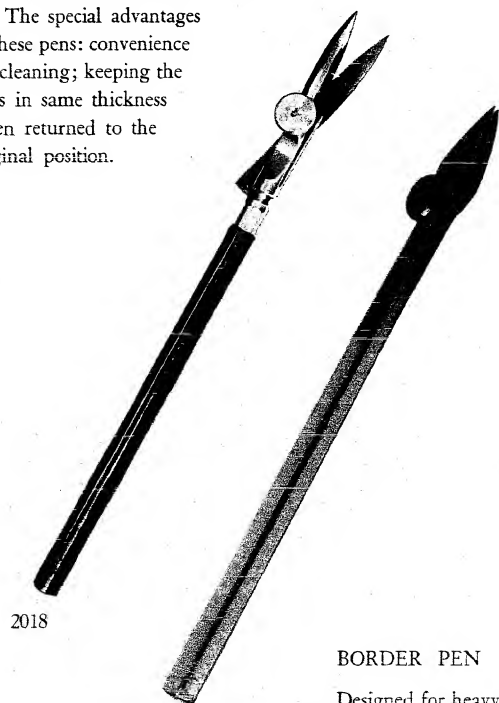2018

1029

BORDER PEN

Designed for heavy and long lines, with aluminium handle.

DOUBLE LINES PEN

Also called "road pen", can be used for drawing very heavy lines when filled with ink between the two pen points.

1035


1033

INDEX PEN

With spring catch and index to the nut to indicate by number the thickness of the line.

CURVE PEN


Specially made for drawing railroads and highways (topographic maps). Not to be used for drawing together with irregular Curve Rule, but for free-hand drawing. It can be also used as ruling pen by turning downward the nut between the handle and the pen.

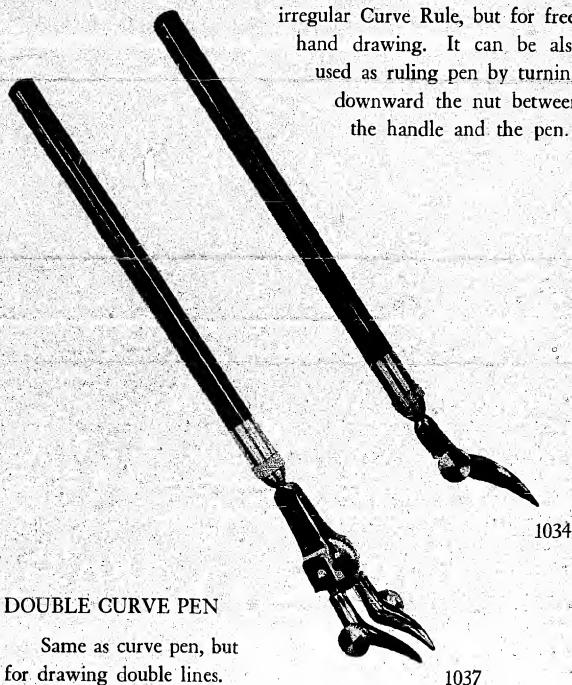DOUBLE CURVE PEN


Same as curve pen, but for drawing double lines.

DROP COMPASS

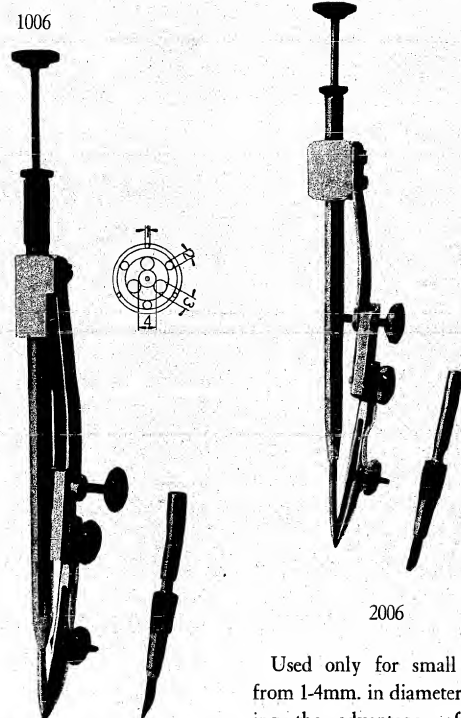Used only for small circles from 1-4mm. in diameter. Having the advantages of time-saving and accuracy.


DOTTING INSTRUMENT

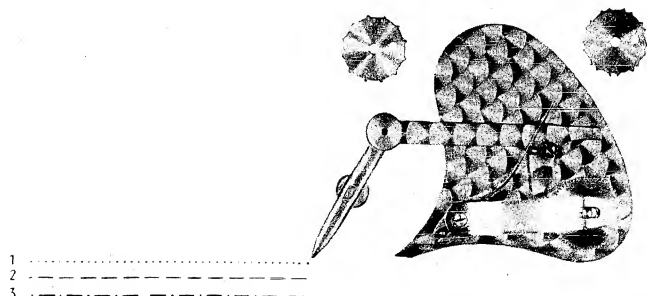1 0 4 2

Specially used for drawing dotting lines, not only saving time but also assuring uniformity. Provided with three dotting wheels for drawing three different kinds of dotting lines. The ruler made in wood of about 3.4mm. thickness is better than that made of any other materials such as metal, plastic and celluloid.

BEAM COMPASS

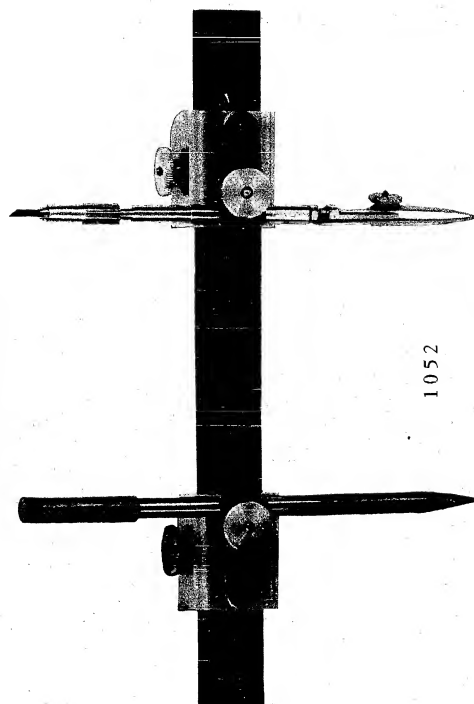
1 0 5 2


For circles of large diameters usually up to 500mm. For special purposes it can be used even up to 4800mm., but in such cases the drawing will require two persons. If the wooden bar is not long enough, the users can replace it with another one of the same thickness and width.

DOTTING INSTRUMENT USED ON BEAM COMPASS

Used for drawing dotting circles by the beam compass, provided with ten dotting wheels.

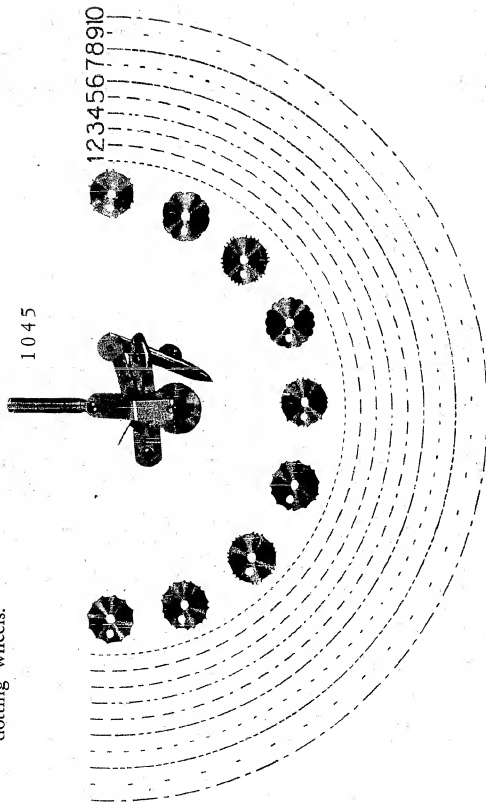PROPORTIONAL DIVIDER

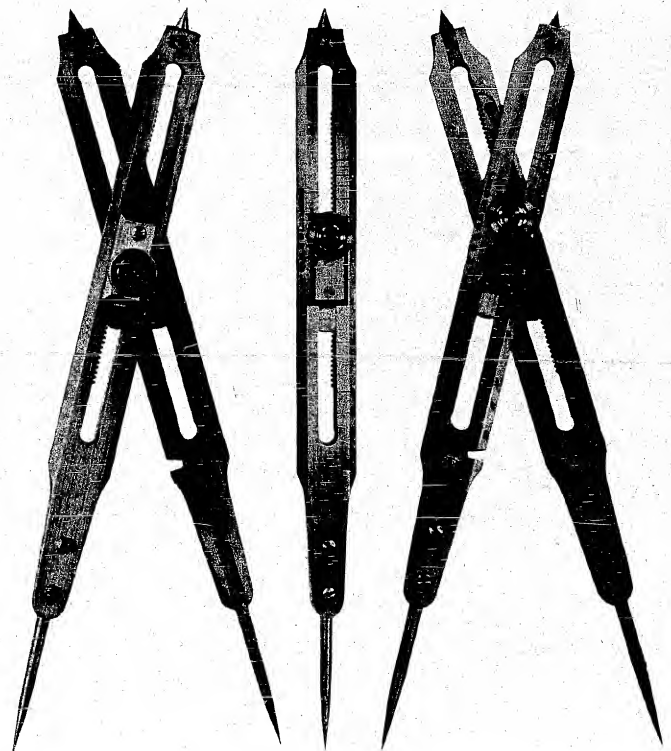1046

7½", with movable points, finely divided for lines, circles, planes and solids, with rack movement.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:


Er-Li-Kou, Hsi Chiao, Peking, China.
Cable Address: "INSTRIMP"

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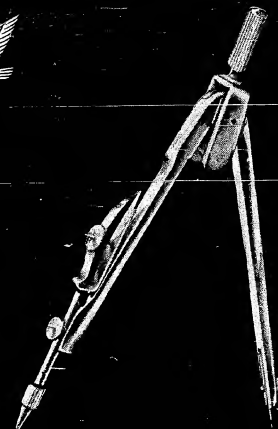
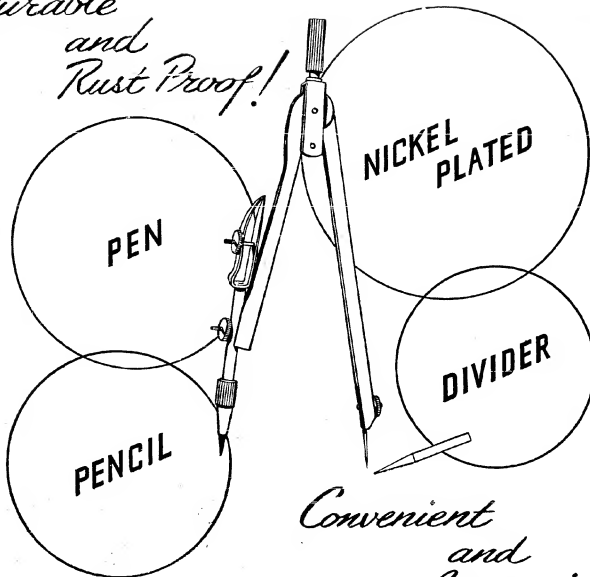
*Convenient
and
Economical!*

SPECIFICATIONS :

OVER ALL LENGTH 110 m.m.
MAXIMUM CIRCLE 110 m.m. in diameter.
PACKING Each in silk lined cardboard box.

THREE PURPOSE COMPASSES

STAT

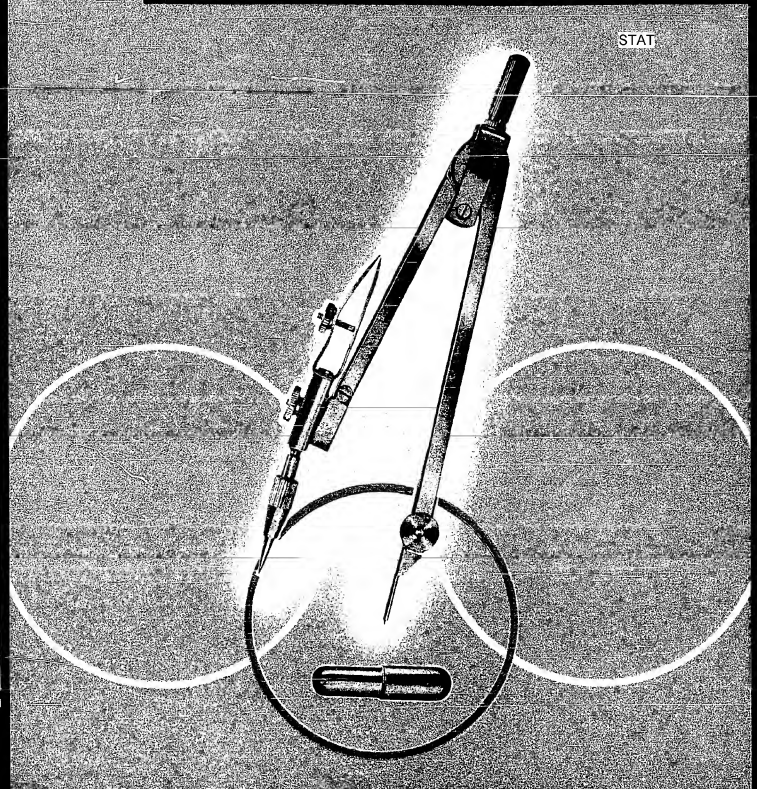CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

THREE PURPOSE COMPASSES


4" compasses with alternative pen, pencil and divider point, suitable for circles below 4" diameter, handsome in appearance and durable in application.

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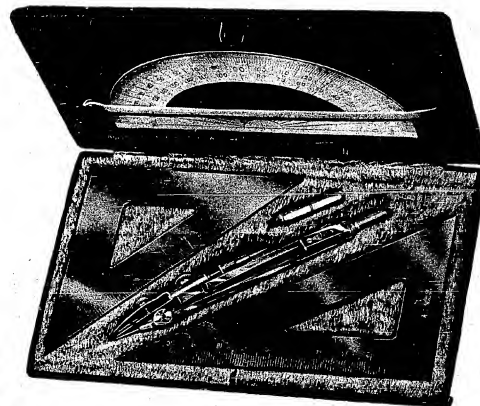CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION
(IMPORTERS & EXPORTERS)

MATHEMATIC INSTRUMENT SET


Attractive students' set consisting of one pair of plastic plain 5" Triangles, one piece of 5" Semi-circular Protractor and one 4" Compasses, with alternative pen, pencil (with lead refill tube) and divider point, suitable for drawing circles below 4" diameter. The whole set being presented in a handsome case with reliable locking device, covered with black Morocco cloth and lined with good quality velvet. It is handy and possesses the advantages of durability and accuracy.

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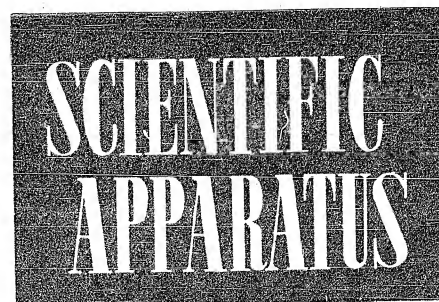Tuning Fork, Vibrograph
Tuning Fork
Tuning Fork
Siren Disc
Articulation Toothed Wheel
Acoustical Resonator
Sonometer

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Diffusion Tester
Brownian Movement Model
Cohesive Force Apparatus
Apparatus showing Transmission of Mechanical Energy into Thermal Energy
Equivalent of Heat
Linear Expansion Indicator
Cubical Expansion Indicator
Double Metal Plates
Heat-conduction Comparator
Indicator for Solid Minification
Indicator for the Maximum Specific Gravity of Water
Radiometer
Indicator for Heat Radiation
Boyle-Mariotte's Law Apparatus
Expansion Apparatus
Illustrator of Capillarity
Crystal space lattice Model
Beam Bending Tester
Testing Apparatus, the Relation Between Temperature and Saturated Atmosphere
Lambrecht's Hygrometer
Dry and Wet Bulb Hygrometer
Hair Hygrometer
Model of Steam Engine, Cross Section Model
Model of Steam Engine, Pumping type
Model of Steam Turbine
Model of Internal combustion Engine
Model of Four-cylinder Internal Combustion Engine
Model of Jet Engine
Air Heat Pressure Coefficient Determination Apparatus

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Glass Rods
Bakelite Rods
Foil Electroscopes
Needle Electroscopes
Electroscope connecting Rod
Electric plume
Electric Streamer
Electric Ball
Coulombmeter
Spherical Electric Conductor

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Induction Cylinder
Cone-shaped Electric Conductor
Metallic Wire Netting Cover
Electrophorus
Influence Machine
Discharge Tongs
Variable Parallel Plate Condenser
Variable Condenser
Leyden jar
Volta Cell
Lead Storage Battery
Resistance Box
Sliding Rheostat
Resistor, Arm Revolving
Constant Value Resistance
Straight Wire Bridge
Galvanometer for Illustration
Suspended Coil Galvanometer
Simple Tangent Galvanometer
Arc Lamp
Testing Apparatus For Current's Chemical Effect
Bar Magnet
Horse-shoe Magnet
Magnetic Needle
Magnetic Particle Model
Iron Filings
Dip Circle
Low Air Pressure Discharge Tube
Low Air Pressure Discharge Tube set
Straight Shot Cathode Ray Tube
Cathode Ray Tube for Magnetic Effect
Cathode Ray Tube for Mechanical Effect
Cathode Ray Tube for Fluorescence
Rontgen Ray Tube
Thermo-electric Couple
Testing device for Mutual Effect of charged conductors
Circular Electro Magnet
Horse-shoe Electric Magnet
Electric Bell
Telegraph Model
Telephone Model
Left hand and Right Hand Rules Indicator
Lenz's Law Apparatus
Self-induction apparatus
Primary and Secondary Coils
Induction Coil
Detachable Transformer

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Universal Transformer
Dynamo Principle Illustrator
Hand driven A.C. and D.C. generator
Motor, 4 Volts
Model of Armature For A.C./D.C. Generator
Illustrator, Magnetic Field Rotary Dynamo
Charging rectifier
Oscillograph
Electric Resonator for Illustration
2-Electrode Electronic Tube
3-Electrode Electronic Tube
Illustrating Plate for Crystal Radio Receiving Set
Illustrating Plate for Single-tube Radio Receiving Set
Universal Electric Measuring Instrument

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Moon and Earth Motion Model
Optical Screen
Full Reflection Illustration Apparatus
Concave Mirror
Double Concave Lens
Double Convex Lens
Lens Set
Prism
Prism, Right-Angle
Micro-Projector
Projector Reflecting Device
Biological Microscope
Microscope Model
Indicating Eye-lens
Dissecting Microscope
Reading Magnifier
Telescope Model
Model of Prism Binoculars
Telescope, Astronomical
Newton's Ring
Color Disc
Colored Glass Plates
Spectroscope
Demonstrating Photoelectric Cell

原子結構學 ATOMIC STRUCTURE

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Wilson's Cloud Chamber
Sphintroscope

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Induction Cylinder
Cone-shaped Electric Conductor
Metallic Wire Netting Cover
Electrophorus
Influence Machine
Discharge Tongs
Variable Parallel Plate Condenser
Variable Condenser
Leyden Jar
Volta Cell
Lead Storage Battery
Resistance Box
Sliding Rheostat
Resistor, Arm Revolving
Constant Value Resistance
Straight Wire Bridge
Galvanometer for Illustration
Suspended Coil Galvanometer
Simple Tangent Galvanometer
Arc Lamp
Testing Apparatus For Current's Chemical Effect
Bar Magnet
Horse-shoe Magnet
Magnetic Needle
Magnetic Particle Model
Iron Filings
Dip Circle
Low Air Pressure Discharge Tube
Low Air Pressure Discharge Tube set
Straight Shot Cathode Ray Tube
Cathode Ray Tube for Magnetic Effect
Cathode Ray Tube for Mechanical Effect
Cathode Ray Tube for Fluorescence
Rontgen Ray Tube
Thermo-electric Couple
Testing device for Mutual Effect of charged conductors
Circular Electro Magnet
Horse-shoe Electric Magnet
Electric Bell
Telegraph Model
Telephone Model
Left hand and Right Hand Rules Indicator
Lenz's Law Apparatus
Self-induction apparatus
Primary and Secondary Coils
Induction Coil
Detachable Transformer

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BIOLOGICAL MODELS, SPECIMENS MICROSCOPIC SLIDES

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VI 浸制解剖標本 PRESERVED AND DISSECTED SPECIMEN

0328	鯽魚解剖	Dissection of Carp (or Silver carp)
0329	青蛙(或蟾蜍)的解剖	Dissection of Frog (or Toad)
0330	龜內臟解剖	Dissection of Turtle
0331	鴿的解剖	Dissection of Pigeon
0332	兔的解剖	Dissection of Rabbit
0333	兔生殖泌尿系統	Excretory and Reproductive System of Rabbit
0334	兔呼吸循環系統	Respiratory and Circulatory organs of Rabbit
0335	脊椎動物五種的腦解剖	Comparison of Brains in the five species of Vertebrates
0336	鯽魚神經系統	Nervous System of silver carp
0337	青蛙(或蟾蜍)的神經系統	Nervous System of Frog (or Toad)
0338	兔的神經系統	Nervous System of Rabbit
0339	鴿的神經系統	Nervous System of Pigeon
0340	蛙的發育順序	Development of Frog
0341	蚌的內部構造	Dissection of clam
0342	蚯蚓解剖(示內臟)	Dissection of Earthworm
0343	蚯蚓解剖(示神經系統和生殖系統)	Nervous and reproductive organs of Earthworm
0345	烏賊解剖示內臟	Dissection of Slug

VII 浸制示範標本 PRESERVED DEMONSTRATION SPECIMENS

0346	海月水母	Jellyfish (<i>Aurelia aurita</i>)
0347	海馬	The sea-horse (<i>Hippocampus</i>)
0348	蛔蟲外形(示雄雌)	<i>Ascaris lumbricoides</i> , male and female
0349	螺	Salamander
0350	蛞蝓	Snail
0351	海星	Starfish
0352	海參	Trepang
0353	比目魚	Flatfish
0354	蜥蜴	Lizard
0355	蜈蚣和馬陸	Scolopendrid and milli pede
0356	狗獾	<i>Esox lucius</i>
0357	蟻	Tapeworm (<i>Taenia</i> sp.)

VIII 剝制標本 MOUNTED SPECIMENS

0358	野鴨(♀)	Wild duck, female
0359	野鴨(♂)	Wild duck, male
0360	啄木鳥	Woodpecker
0361	蒼鷹	Hawk
0362	貓頭鷹	Owl
0363	鴿	Pigeon
0364	北京鴨	Peking duck
0365	蝙蝠	Bat
0366	松鼠	Squirrel
0367	鼹鼠	Mole

IX 昆蟲標本 INSECTS

0369	蝗蟲生活史	Life history of Grasshopper (<i>Locusta migratoria</i>)
0370	菜粉蝶生活史	Life history of the Cabbage butterfly (<i>Pieris rapae</i> L.)
0371	柞蠶生活史	Life history of the Chinese oak silkworm (<i>Antheraea Pernyi</i> Guesin)
0372	蜜蜂生活史	Life history of Honey-bee (<i>Apis mellifera</i> L.)
0373	蠶蛾生活史	Life history of Silkworm (<i>Bombyx mori</i> L.)
0374	蠅的生活史	Life history of Fly (<i>Musca domestica</i>)
0375	螟蛾生活史	Life history of Corn borer (<i>Pyrausta nubilalis</i> Hubner)
0375	竹節蟲	Walking-stick (<i>Diapheromera orientalis</i>)
0376	枯葉蝶	Leaf butterfly

X 顯微鏡玻片標本 MICROSCOPIC SLIDE SETS

BZ 25	動物學	Zoology	25	pcs.	set
BZ 50	"	"	50	"	"
GB 25	植物學	Botany	25	"	"
GB 50	"	"	50	"	"
GH 25	生理組織學	Physiology & Histology	25	"	"
GH 50	"	"	50	"	"
PH 12	寄生蟲學	Parasitology	12	"	"
PH 25	"	"	25	"	"
BS 25	生物學	Biology	25	"	"
BS 50	"	"	50	"	"
BM 25	初中示範教	Assorted set for Junior Middle School	25	"	"
BM 50	"	"	50	"	"

Exported by:

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

HEAD OFFICE:

En-Li-Kou, Hsi Chiao, Peking, China
Cable Address: "INSTRIMPOT"

BRANCH OFFICE:

SHANGHAI: 27 Chungshan Road (E. 1)
Cable Address: "INSTRIMP"

TIENSIN: 171 Kien Shieh Road,
Cable Address: "INSTRIMP"

CANTON: 25 Tai Ping Road, S.
Cable Address: "INSTRIMP"

AGENCY:

HONGKONG: CHINA RESOURCES CO.
Bank of China Building 11/12 Fl.
Cable Address: "CIRECO HONGKONG"

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION^{STAT}

EXPORT LIST

1. Radio Equipments:

RADIO RECEIVER: 5-tube super-heterodyne radio receiver with medium and short wave bands in beautiful bakelite and/or veneer wood cabinet, having the best features of high sensitivity and selectivity giving a perfect tonal reproduction. Power supply: 110/220 volts A.C. 50 cycles.

Radio Accessories:

P.M. Loud Speakers: 5", 6½" and 8".

Radio Tubes: 6SA7GT, 6SK7GT, 6SQ7GT, 6V6GT, 5Y3GT.

I.F. Transformer: Air Core. Iron Core.

Radio Coils: Antenna, OSC., R.F. Long and short wave coils.

Tubular Paper Condensers: 0.0001, 0.00025, 0.0005, 0.001, 0.002, 0.005, 0.006, 0.01, 0.02, 0.05, 0.1, mfd.

Dry Electrolytic Condensers: 8, 25, 30, mfd.

Variable Condensers: 2-gang, 3-gang. .00035 mfd.

Compound Radio Batteries: A-battery 1.5 v. 39 ohm.
B-battery 91.2 v. 4800 ohm.

Toggle Switches: S.P.S.T. S.P.D.T. D.P.S.T. D.P.D.T. and Push button switch.

Head Phone: 3000 ohm.

Other Parts: Phone plug and jack, 8-prong sub-panel socket, Crystal unit and detector, Padding condenser.

Soldering Irons: 20-45-75-100-150. W. 110 V. or 220 V.

Resistors: Wire wound or carbon.

Phono Equipments:

Automatic record changer with magnetic pickup, can be supplied with or without Amplifier, 110 or 220 V. AC. 50 cycles.

Electric Record Player with magnetic pickup, 110 or 220 V. AC. 50 cycles.

Permanent Magnet pickup cartridges.

Public Speech Power Amplifiers: 15 W. 25-30 W. with or without medium & short wave receiving equipments.

Microphones: High impedance with the dynamic or ribbon.

Tape and Wire Recorders.

2. Photographic Equipments for Photo Studio:

Studio Portrait Camera,	No. 8 for	6 1/2" x 8 1/2"	Negatives
Complete Outfit:	No. 12 "	10" x 12"	"
Portable View Camera,	No. 4 "	3 1/2" x 4 1/2"	"
Complete Outfit:	No. 6 "	4 1/2" x 6 1/2"	"
	No. 8 "	6 1/2" x 8 1/2"	"
	No. 10 "	8" x 10"	"
	No. 12 "	10" x 12"	"

Photographic Enlarger: No. 4 for negatives up to 3 1/2" x 4 1/2"
No. 6 for negatives up to 4 1/2" x 6 1/2"

Retouching Cabinet: No. 8 for 6 1/2" x 8 1/2" negatives
No. 12 for 10" x 12" negatives

Spot Lights: 6" or 8" Condenser, 750W-1000W.
with or without Revolving Disc.

Glossy Print Dryer, single or double side. Size: 11" x 14"
Power Consumption: 1000W, 110 or 220V. 14" x 18"
18" x 24"

Photo-Flood Light-units: 1, 2, 3, 4, 5, 6-light units of different types.

Universal Tripod: Overall dimensions 46"

Range of elevation 20"

Universal Adapter: No. 4 for Negative of 3 1/2" x 4 1/2"

No. 6 for Negative of 4 1/2" x 6 1/2"

Special kits are provided for. Capable of taking pictures of different sizes at different intervals.

Universal Film Adapter: Specially made for No. 120 or 620 Films.
Capable of taking 11 Exposures. Kits are also provided for exposures of half and quarter sizes.

Printing Frames, all sizes.

Negative Drying Racks, all sizes.

Contact Printers, all sizes.

Negative Holders, all sizes.

Photographic Contact Printing & Enlarging Paper.

Photographic Carbon: 13 x 305 mm.

3. Electrical Measuring Instruments:

Precision AC/DC Ammeters and Voltmeters, portable.

Accuracy $\pm 0.5\%$, 12 Ranges, Dimension 264 x 216 x 158 mm.

Precision Single and Three Phase Wattmeter, portable.

Accuracy $\pm 1\%$, 24 Ranges, Dimension 264 x 261 x 147 mm.

Precision Single and Three Phase Power Factor Meter, portable.

Accuracy Cos. 1 ϕ , 9 ranges, Dimension 201 x 261 x 147 mm.

Precision Cycle Counter, portable.

Rated Voltage: 110-220 Volts, Dimension 175 x 155 x 100 mm.

Precision Resistance Tester (Megohmmeter), portable.

Measuring ranges 0-20-200 Megohms, Dimension 255 x 210 x 205 mm.

Precision Kelvin Bridge; Wheatstone Bridge; Potentiometer.

Capacitance Box; Resistance Box; portable.

Switchboard Voltmeter, Ammeter; AC/DC 4", 7" and other sizes.

Switchboard Power Factor Meter, Wattmeter.

Household Current Meter.

4. Cinematographic Equipments:

16 mm. Sound Film Projector. Complete with accessories.
Slide Projectors.

Cinema Carbons: H.I. or L.I. 8 x 305 mm, 7 x 230 mm. 7 x 305 mm.
6 x 230 mm. Other sizes made to order.

Projection Lamps.

Exciter Bulbs.

5. Telephone Equipments:

Desk and Wall Telephones.

Automatic Telephone Equipments.

Telephone Equipments, Common Battery System.

Switchboards.

Alarm Receiving & Despatching Apparatus.

6. Scientific Instruments:

Yarn Tensile Strength Tester for Single Yarn.

Cotton Drawing Apparatus } For Cotton Maturity, Length and
Cotton Sorter } Fineness Tests.

Motor Driven Lap Meter.

Motor Driven Wrap Reel.

Cloth Strength Tester (Hydraulic Type).

Lea Yarn Strength Tester (Hydraulic Type).

Roving & Sliver Tester (Motor Driven).

Shirley Analyzer.

Twist Tester.

CO₂ Analyzer with Remote Control and Multi-Indication Meter.

Flowmeter, Electrical and Mechanical, with Remote Control and
Automatic Recording Devices, for Gases and Liquids.

Water Meter.

Planimeter: For Calculation of Areas of Plane Figures.

Pyrometers: Thermo-Couple and Optical.

Metallurgical Microscope.

Micro Projection Equipments.

7. Laboratory Equipments:

A. Analytical and Precision Balances:

Analytical Balances:

1. Single Pan Projection Reading:
Capacity: 100 gms.
Accuracy of Weighing: by difference
within the Optical Range: 0.05 mg.

2. Automatic Projection Reading with
Air-Damping Device for rapid and reliable weighing.

Capacity: 200 gms.
Sensitivity: 0.1 mg.

Precision Balances:

for Schools and Laboratory Experiments.
Sensitivity from 0.4 to 20 mg.

- B. Glass and Porcelain Laboratory Apparatus.
- C. Chemical Thermometers, all kinds.
- D. Magnifiers: 2x to 20x.
- E. Hydrometers.
- F. Spirit Levels.

8. Educational Equipments:

Physics Appliances for Schools.

Anatomical and Biological Models.

Biological Microscopes: 675x, 1500x.

Biological Dissecting Microscopes: 20x.

Biological Microscopic Slides.

Precision Drawing Instruments: 4 - 28 pcs. per set.

3-Purpose Compass.

Mathematical Instrument set,

Comprising, of 1 3-purpose

Compass, 1 Plastic protractor

and 1 pair of plastic

Triangles.

For Junior Schools.

Terrestrial Globes.

9. Office Appliances:

Calculating Machines.

Chinese Typewriter (Universal) and Accessories.

Household Thermometers. Wall, Dry and Wet, Max. and Min.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

27 Chungshan Road (E. 1),
Shanghai China.

Cable Address: "INSTRIMP"

(Head Office at Peking, Branch Offices
at Tientsin & Canton)

Biological Microscopic Slides.

Precision Drawing Instruments: 4 - 28 pcs. per set.

3-Purpose Compass.

Mathematical Instrument set,
Comprising, of 1 3-purpose
Compass, 1 Plastic protractor
and 1 pair of plastic
Triangles.

For Junior Schools.

Terrestrial Globes.

9. Office Appliances:

Calculating Machines.

Chinese Typewriter (Universal) and Accessories.

Household Thermometers, Wall, Dry and Wet, Max. and Min.

CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

27 Chungshan Road (E. 1),

Shanghai China.

Cable Address: "INSTRIMP"


(Head Office at Peking, Branch Offices
at Tientsin & Canton)

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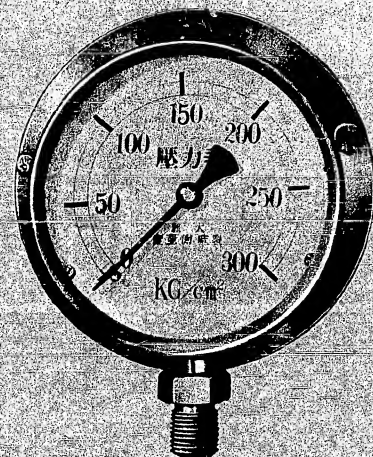
CHINA NATIONAL INSTRUMENTS IMPORT CORPORATION

GENERAL USE: Precision pressure measuring apparatus for steam boilers of all sizes, vehicle oil compressors, air compressors etc. The measurement result can be easily read by the scales indicated.

APPLICATIONS: A. Suitable for measuring pressure of slow flow and even liquid and atmosphere, but not suitable for measuring explosive matters and crystallized substances.

B. The temperature of the liquid and atmosphere to be measured should range from -20°C to 60°C .

PARTCUAPS:

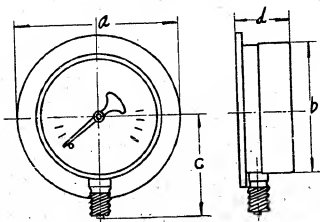(unit=mm)

fixed specifications (Diameter)	a (dia)	b (dia)	c	d
60	—	60	51	36
75	98	79	59	41
100	126	102	76	41
125	155	128	127	49
150	180	154	127	48
250	—	256	161	54

MAIN SPECIFICATIONS


specification of outside circumference (unit:mm)				minimum degreeind- icated	Distance between two indicated degrees	Permissible pressure	Permis- sable tolerance (±)
60	75~100	125~150	250				
Pressure limit							
Unit: Kg/cm ²							
0.5	0.5	—	—	0.02	0.1	0.34	1/2 of the smallest indicated degree
1	1	—	—	0.05	0.2	0.7	
2	2	2	—	0.05	0.5	1.35	
3	3	3	—	0.1	0.5	2	
4	4	4	—	0.1	1	2.7	
5	5	5	—	0.2	1	3.4	
6	6	6	—	0.2	1	4	
7	7	7	—	0.2	1	4.8	
8	8	8	—	0.2	2	5.4	
10	10	10	10	0.5	2	7	
12	12	12	12	0.5	2	8	
16	16	16	16	0.5	2	11	
20	20	20	20	0.5	5	13.5	
25	25	25	25	1	5	17	
30	30	30	30	1	5	20	
35	35	35	35	1	5	24	
40	40	40	40	1	10	27	
50	50	50	50	2	10	34	
60	60	60	60	2	10	40	
80	80	80	80	2	10	54	
100	100	100	100	5	20	70	
120	120	120	120	5	20	80	
160	160	160	160	5	20	110	
200	200	200	200	5	50	135	
250	250	250	250	10	50	170	
—	300	300	300	10	50	200	
—	—	400	400	10	50	270	
—	—	500	500	20	100	340	
—	—	600	600	20	100	400	
—	—	800	800	20	100	540	
—	—	1000	1000	50	200	700	
—	—	—	—	—	—	—	

INSTALLATION:


- A. If the gauge is to be installed near the boiler with a surface temperature within 100°C . It should be equipped with curved tubes.
- B. Should the surface temperature be over 100°C , tubes of proportionate length is to be installed.
- C. The distance from the surface of the gauge to the boiler should be 250-300mm.

POSITION OF INSTALLATION

Picture 1

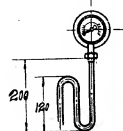Picture 2

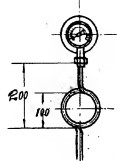- REMARKS: A. Performance under even pressure loading, pressure is allowed to the extent to $2/3$ of the pressure limit.
- B. Performance under unsteady pressure loading, pressure is allowed to the extent of $1/2$ of the pressure limit.
- C. If the performance of air pressure gauge is to be over 100 kg/cm^2 , oils should not be emerged into the tube, but this is not applicable to oil compressors.

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CHINA NATIONAL INDUSTRIAL

PRODUCTS EXPORT CORPORATION TIENTSIN BRANCH

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Flasks, wide mouth, with flat bottom, short neck

Capacity: ml.

100 150 250 500

Price: @SH.

0/11 1/3 1/6 1/10

CHINA NATIONAL INSTRUMENTS IMPORT CORP.
SHANGHAI
(IMPORTERS & EXPORTERS)

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CHEMICAL APPARATUSPRICE LIST

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Flasks, with round bottom

Capacity: ml.	50	100	150	200	250	300	500	1000	1500	2000
Price: @SH.	0/8	0/9	0/11	1/1	1/3	1/5	1/8	2/7	4/7	6/5
	2500	3000	5000	10000	15000	20000	25000			
	7/11	9/10	15/2	28/9	42/5	54/6	63/8			

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Flasks, with round bottom & rolled neck

Capacity: ml.	250	500	1000	1500	2000	2500	3000	5000	7500
Price: @ SH.	1/6	2/-	3/-	5/2	6/5	7/11	9/10	15/2	21/-
						10000	15000	20000	25000
						28/9	42/5	54/6	63/8

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Flasks, with round bottom & short rolled neck

Capacity: ml.	250	500	1000	1500	2000	2500	3000	5000	7500
Price: @ SH.	1/6	2/-	3/-	5/2	6/5	7/11	9/10	15/2	21/-
						10000	15000	20000	25000
						28/9	42/5	54/6	63/8

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Retorts with tubulature and ground in stopper

Capacity: ml.	125	250	500	1000	2500	5000	7500	10000
Price: @ SH.	4/7	6/-	7/9	9/10	19/8	33/4	48/6	68/2

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Length of tube: mm. 150 180 200

Price: @ SH. 0/8 0/10 1/2

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Connection tubes with branch, T Shaped

Price: @ SH. 0/3

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Connection tubes with branch, Y Shaped

Price: @ SH. 0/4

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Connection tubes with branch, U Shaped

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Aspirator bottles ground in stopcocks & stopper

Capacity: ml.	2500	5000	10000	20000
Price: @ SH.	18/2	25/9	37/11	72/8

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Funnel, short stem

Diameter: mm.	40	50	60	75	90	120	150	180	210	240	300
Length of tube mm.	50	50	60	75	90	120	150	150	180	180	180
Price: @ SH.	0/9	0/11	1/1	1/3	1/5	1/8	2/2	3/-	4/7	7/7	12/1

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Funnel, angle 60°, long stem

Diameter: mm.	50	60	75
Length:	150	150	150
Price: @ SH.	1/6	1/10	2/3

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Funnel tubes, thistle top, straight form

Length of tubes: 300 mm. Price: @ SH. 0/7

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Funnel tubes, thistle top, with loop

Length of tubes: 300 mm. Price: @ SH. 0/9

POOR ORIGINAL

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CHEMICAL APPARATUSPRICE LIST

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Funnel tubes, thistle top with loop and bulbs

Length: 500 mm.

Price: @ SH. 0/11

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Funnel tubes, thistle top with double loop and 2 bulbs

Length: 300 mm.

Price: @ SH. 1/3

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Separating funnels heavy globe shape

Capacity: ml. 1000 2000 3000 5000

Price: @ SH. 30/4 45/5 68/2 90/11

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Graduated Cylinders, with spout

Capacity: ml. 10 25 50 100 250 500 1000 2000

Graduation: ml. 1/5 1/2 1 2 5 5 10 20

Price: @ SH. 1/6 1/10 2/2 2/7 4/7 6/10 9/10 18/-

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Stoppered graduated cylinder

Capacity: ml. 10 25 50 100 250 500 1000 2000

Graduation: ml. 1/5 1/2 1 2 5 5 10 20

Price: @ SH. 2/6 3/4 3/11 4/10 6/10 9/10 13/8 21/-

POOR ORIGINAL

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CHEMICAL APPARATUSPRICE LIST

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Graduated cylinder, for acids

Price: @ SH.

1/6

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Graduated cylinder conical

Capacity: ml. 10 20 50 100 250 500 1000

Graduation: ml. 1 2 5 5 10 20 50

Price: @ SH. 1/6 1/10 2/2 2/7 4/7 7/7 10/7

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Volumetric flasks with stopper

Capacity: ml. 25 50 100 250 500 1000 2000

Price: @ SH. 3/8 3/11 4/7 5/4 6/- 9/1 13/8